

# Illumina Sequencing Overview: Library Prep to Data Analysis

Ryan Gentry M.S.  
Sequencing Specialist



QB7845

**For Research Use Only. Not for use in diagnostic procedures.**

© 2019 Illumina, Inc. All rights reserved.

Illumina, 24sure, BaseSpace, BeadArray, BlueFish, BlueFuse, BlueGnome, cBot, CSPro, CytoChip, DesignStudio, Epicentre, ForenSeq, Genetic Energy, GenomeStudio, GoldenGate, HiScan, HiSeq, HiSeq X, Infinium, iScan, iSelect, MiniSeq, MiSeq, MiSeqDx, MiSeq FGx, NeoPrep, NextBio, Nextera, NextSeq, Powered by Illumina, SureMDA, TruGenome, TruSeq, TruSight, Understand Your Genome, UYG, VeraCode, verifi, VeriSeq, the pumpkin orange color, and the streaming bases design are trademarks of Illumina, Inc. and/or its affiliate(s) in the US and/or other countries. All other names, logos, and other trademarks are the property of their respective owners.

**illumina®**

# Comparing Technologies



## q/RT-PCR

q/RT-PCR allows for the analysis of particular variants at specific locations

### + Benefits

- High sensitivity
- Capital equipment already found in most labs

### - Limitations

- Low discovery power
- Low variant resolution
- Low scalability

## Sanger/CE

Sanger/CE is able to interrogate a gene of interest

### + Benefits

- Cost effective for small stretches of DNA
- Well known technique

### - Limitations

- Low sensitivity (down to 20%)
- Low discovery power
- Low scalability

## Targeted NGS

Targeted NGS allows for simultaneous screening of several hundreds to thousands of genes

### + Benefits

- Expanded discovery power
- Maintain resolution and high sensitivity
- More data from smaller amounts of DNA and RNA
- Higher throughput with sample multiplexing


### - Limitations

- May be less cost effective when interrogating a low number of samples

# Illumina Sequencing Workflow

1

## Library Preparation


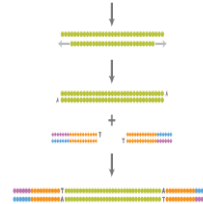


3

## Sequencing

[illegible]

# Library Prep is Critical for Successful Sequencing

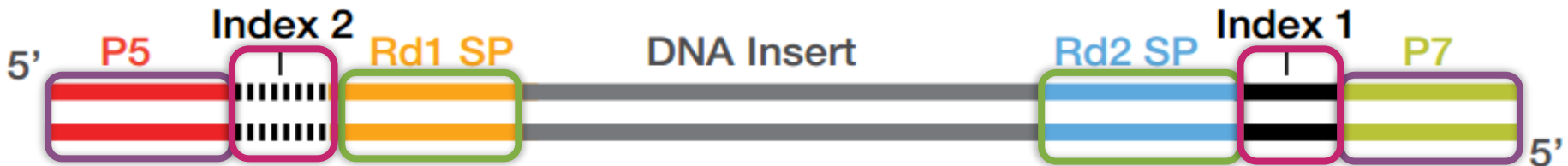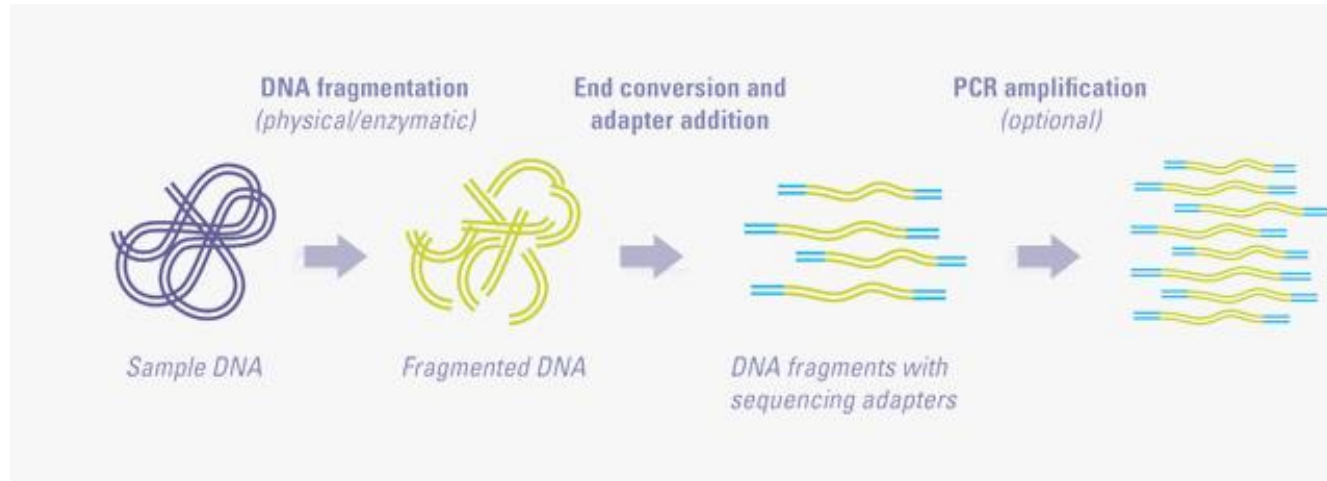

**For clustering:**  
*Libraries must have P5 and P7 binding regions on either end of a library*

**For sequencing:**  
*Libraries must have sequencing primer binding regions*

**For mixing samples:**  
*Libraries must have a unique index or barcodes sequence*


# AmpliSeq™ for Illumina® workflow



Panel



Index plate

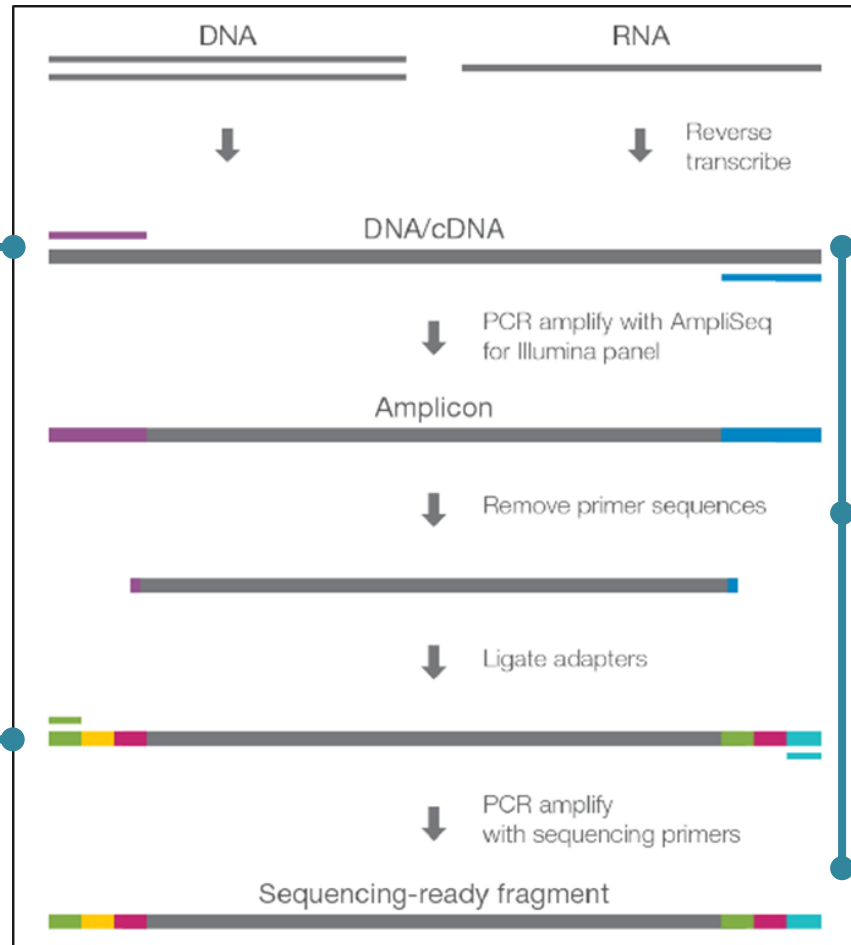

Library PLUS


1–100 ng Input  
Recommend – 10ng / Pool

5 – 7 hrs Total Assay Time


1.5 hrs Hands on Time

# Tools for DNA Library Preparation

*Fast or high performance*

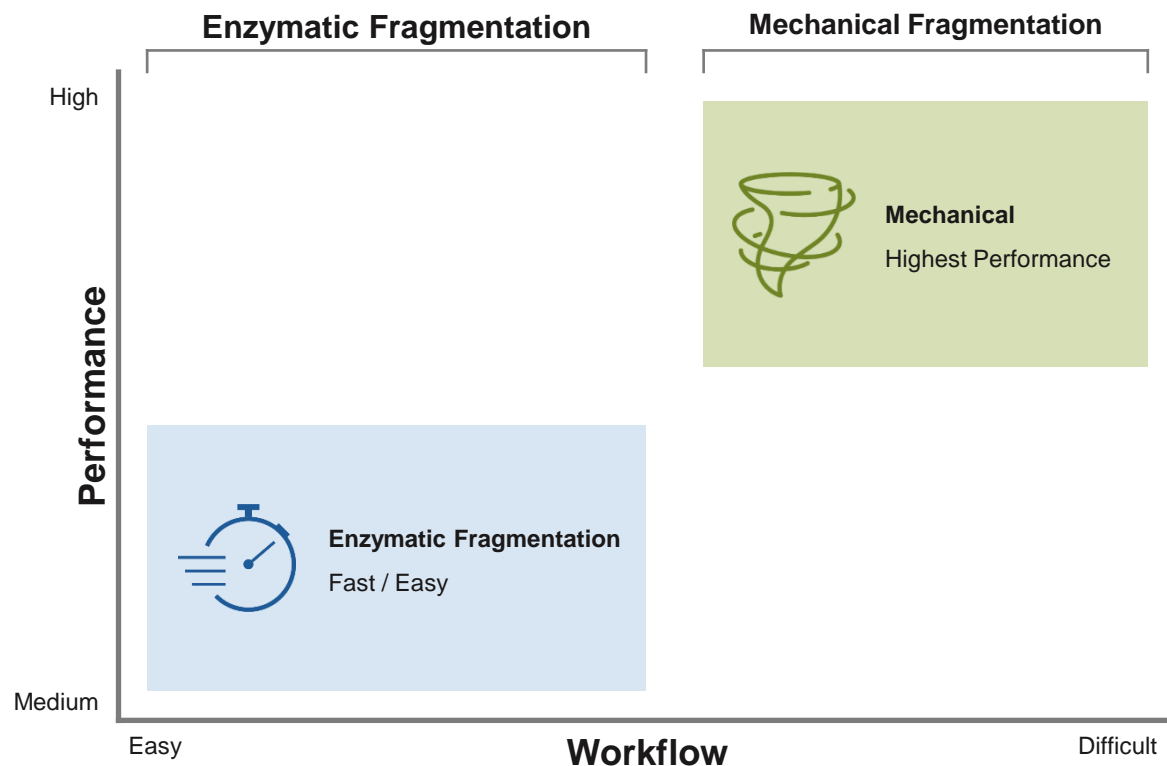

# Mechanical fragmentation workflow

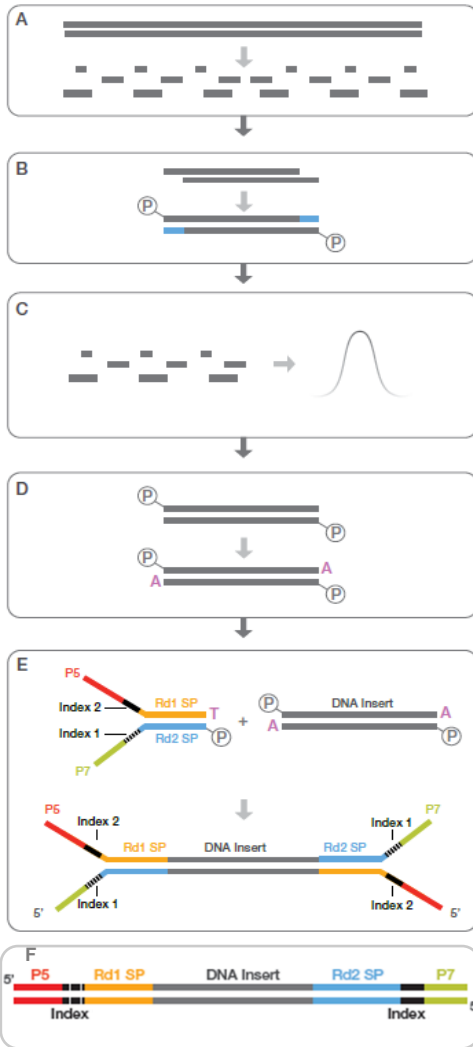

A. Genomic DNA is fragmented

B. Genomic DNA is fragmented


C. Fragments narrowly size selected

D. A-base added

E. Adapters ligated

F. (only for TruSeq Nano kit) PCR enriches complete libraries


# Enzymatic fragmentation workflow

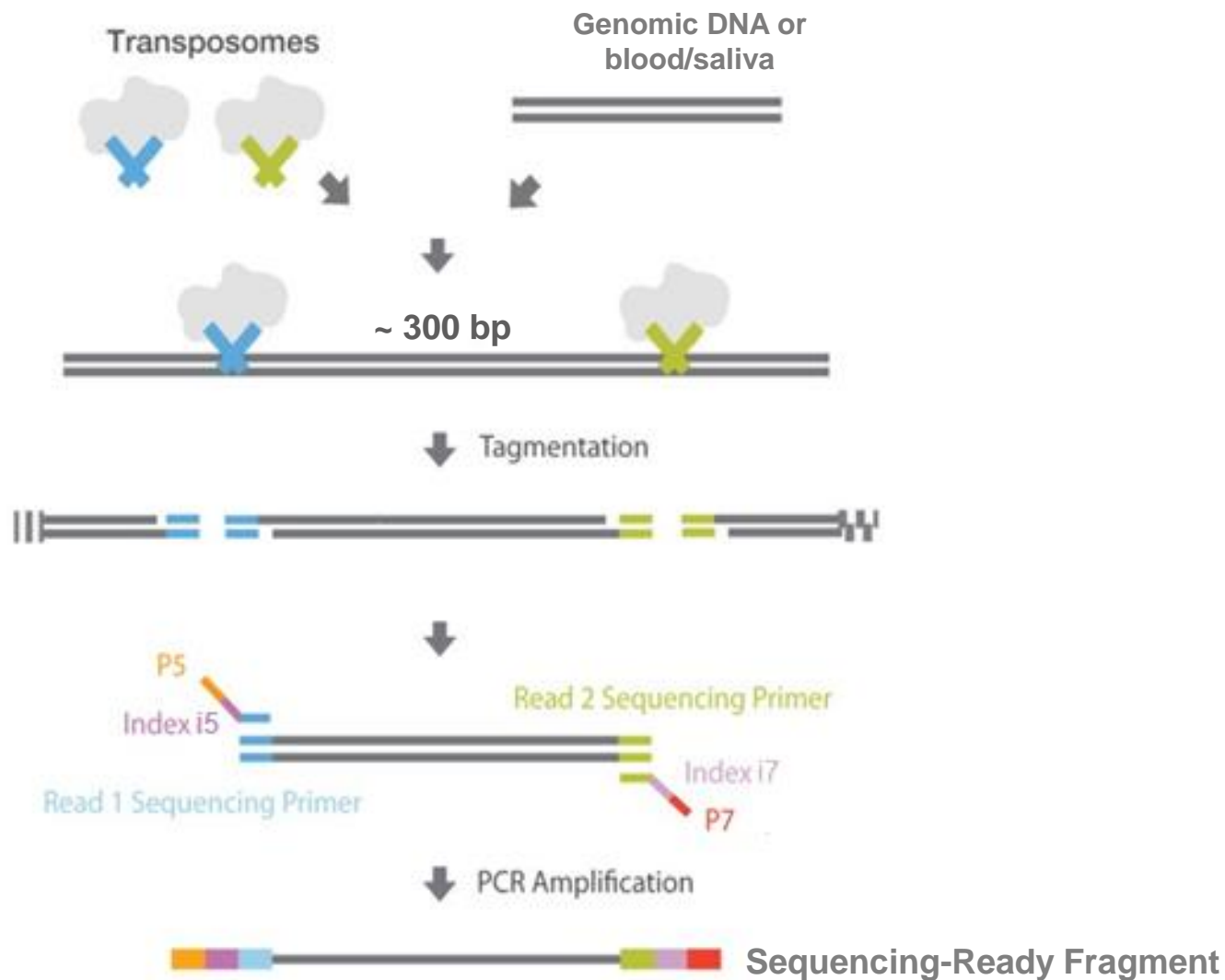




# Introducing Illumina DNA Prep, (M) Tagmentation

*One DNA prep, multiple solutions*

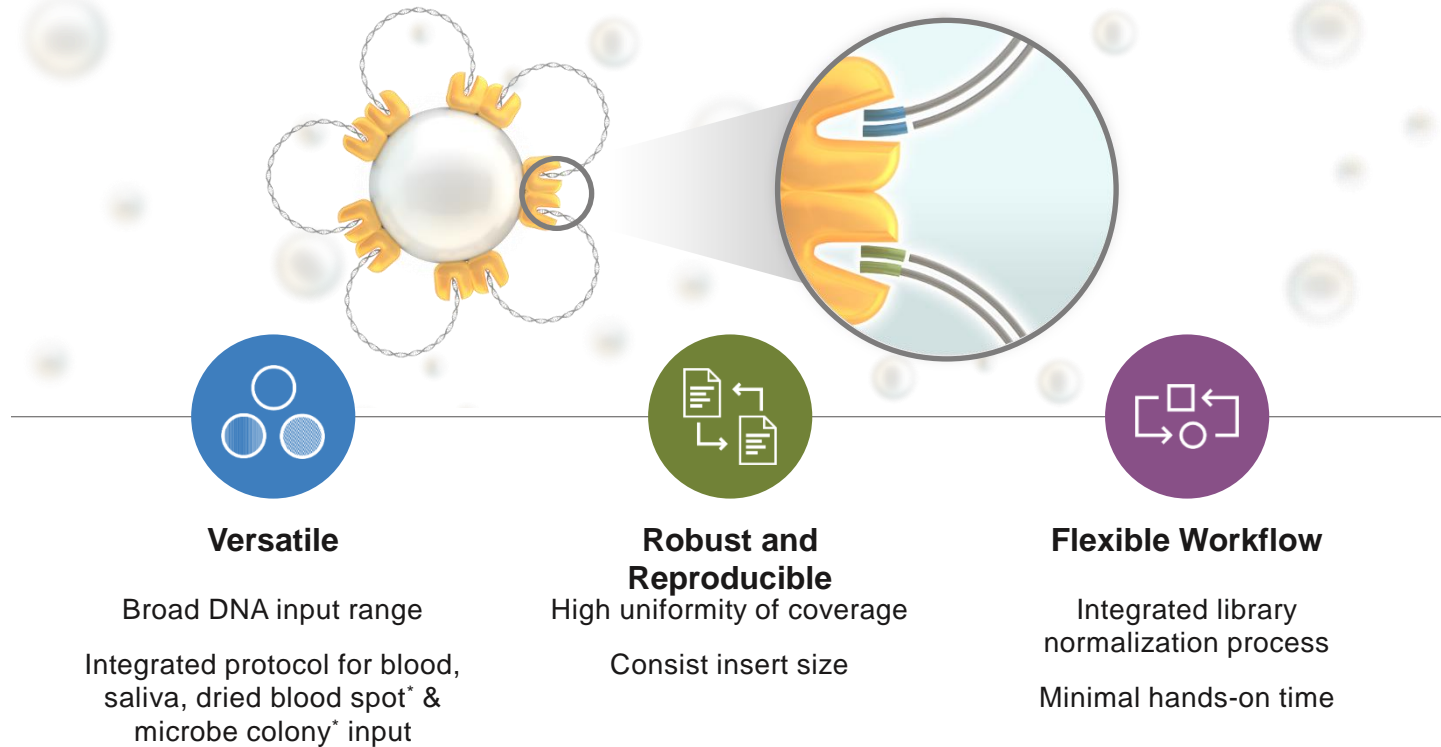


\*Demonstrated protocol

# Rapid Capture Enrichment Workflow

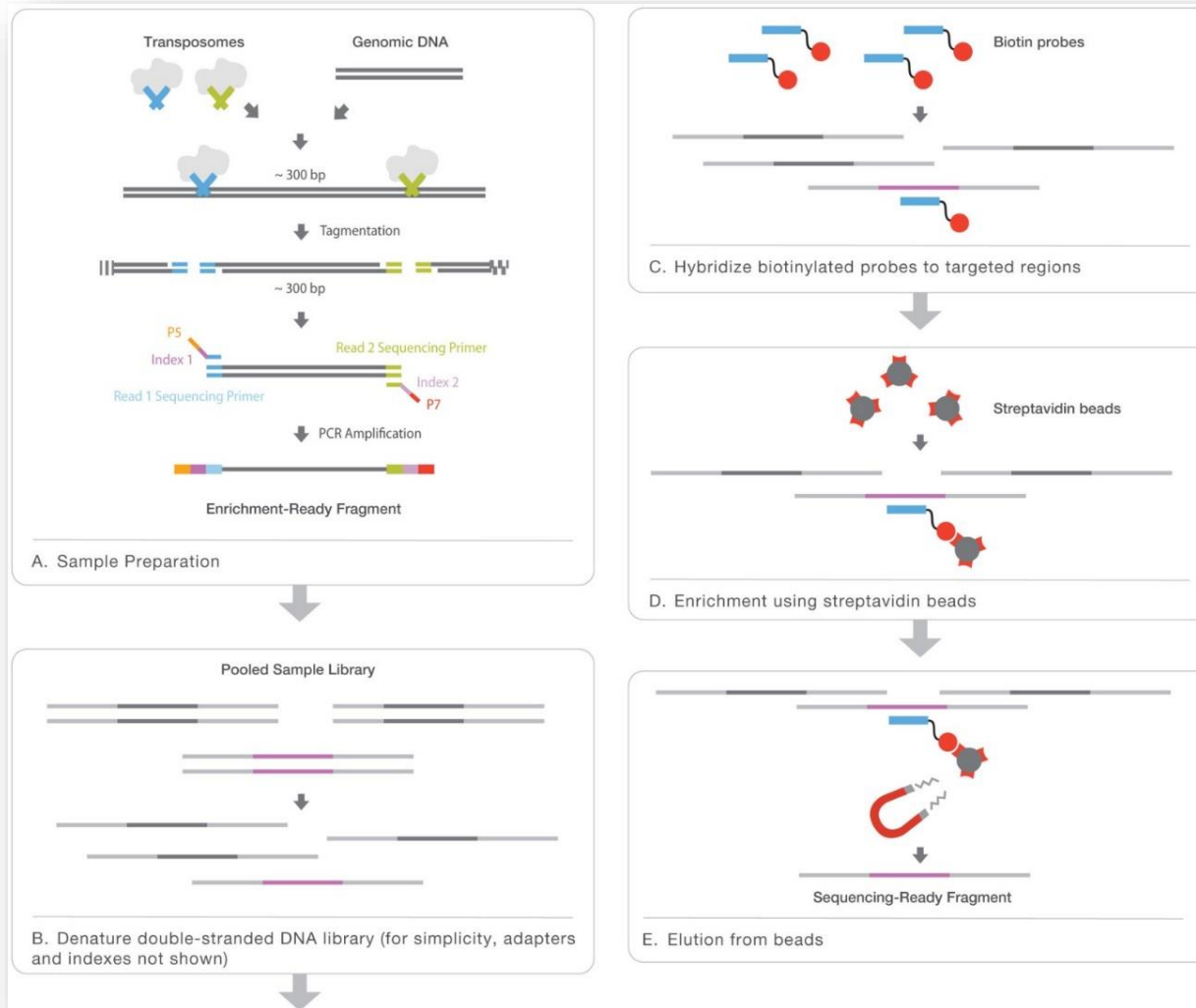


# One Technology, Many Applications

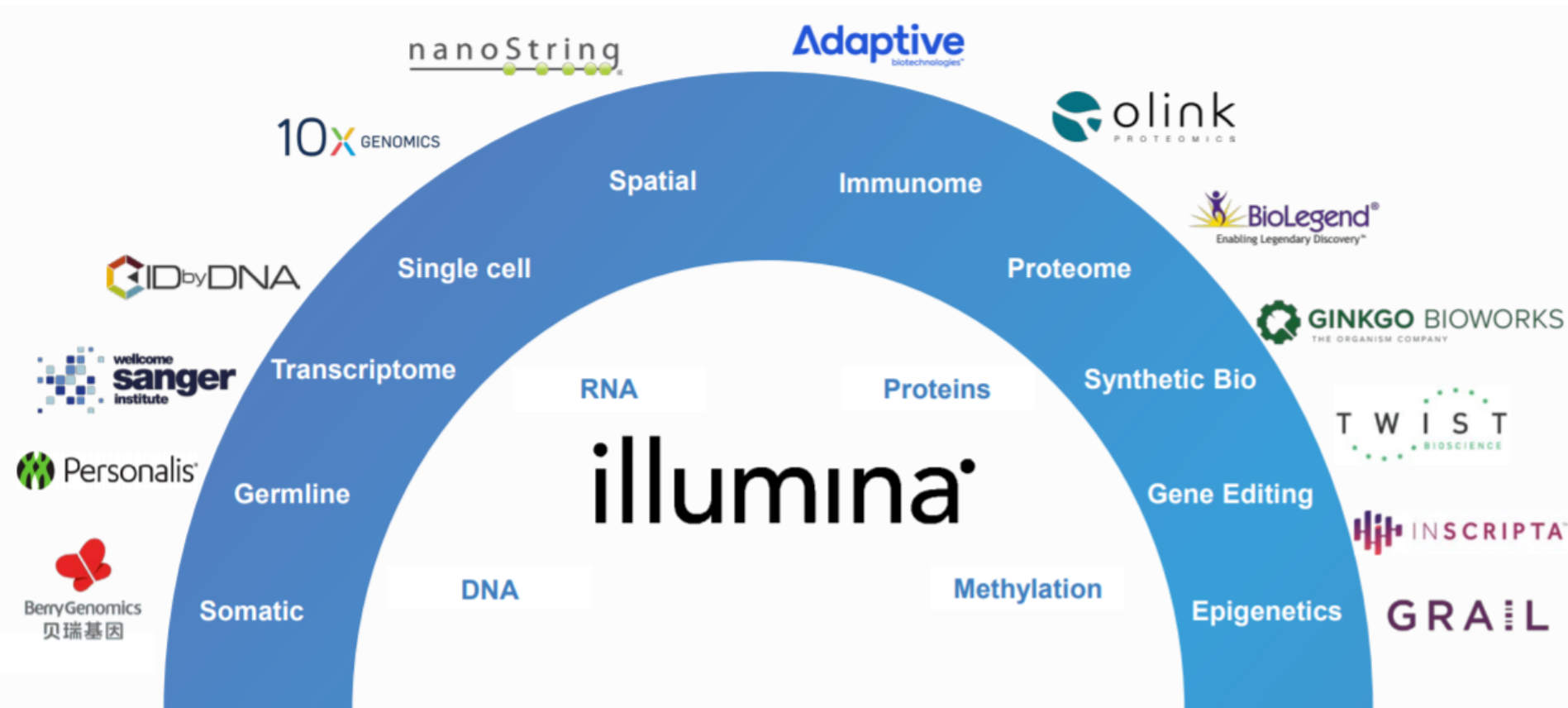

# Illumina Sequencing Workflow


1

Library Preparation

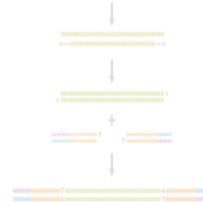


2

Cluster Generation

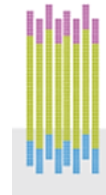

3

Sequencing

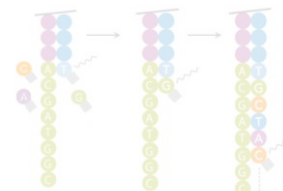


4

Data Analysis

```
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
5TAAGGCTAGGTTTCATGCTA
```

# Illumina Sequencing Portfolio

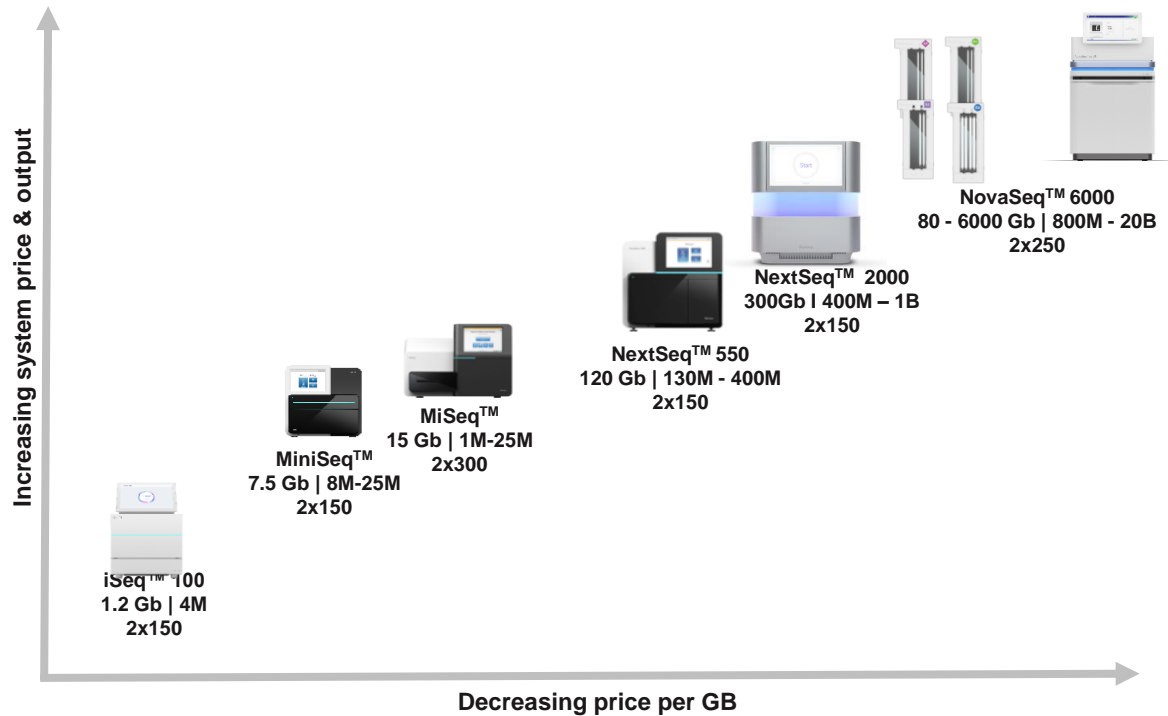

For Research Use Only. Not for use in diagnostic procedures.

# MiSeq System



PE300 | **75% > Q30**

---

**0.3–15** | gigabases

---

**1–25 million** | clusters

# NextSeq 550 has Tunable Output, High Data Quality, and Array Capabilities



## High-Output

Up to 120 Gb  
400M clusters PF  
1x75 bp to 2x150 bp



## Mid-Output

Up to 40 Gb  
130M clusters PF  
2x75 bp to 2x150 bp



## MethylationEPIC and CytoSNP 850K

Clinically relevant content on arrays



## SBS chemistry

90% of all NGS data are generated on an Illumina platform  
High accuracy  
High quality



## DRAGEN™ pipelines on BaseSpace™

Industry-leading accuracy and speed  
Variants include indels, small variants, and CNVs  
In the cloud with BaseSpace and on prem with a server

For Research Use Only. Not for use in diagnostic procedures.

# MiSeq Offers Scalable Sequencing



1 Million Reads

4 Million Reads

15 Million Reads

25 Million Reads

## MiSeq

### Core Consumables Version 2 Nano

- 500 cycles (Nano)
- 300 cycles (Nano)

### Core Consumables Version 2 Micro

- 300 cycles (Micro)

### Core Consumables Version 2

- 500 cycles
- 300 cycles
- 150 cycles

### Core Consumables Version 3

- 600 cycles
- 150 cycles



# NextSeq 1000/2000 Configurations

## NextSeq 1000

\$210K | Now Available



### NextSeq 1000

120GB

Max output  
(100/200/300 cycles)

### Field

**Upgradeable**

To NextSeq 2000  
\$150K

## NextSeq 2000

\$335K | Now Available



### P2 Flow Cell

120GB

Max output  
(100/200/300 cycles)

### P3 Flow Cell

**330GB**

Max output  
(50/100/200/300 cycles)

## Further workflow improvements

Custom primers and  
custom recipe

DRAGEN Single Cell RNA,  
DRAGEN Enrichment Somatic


FastQ compression and  
FastQC metrics

Additional workflow features /  
capabilities and BSSH / control  
software improvements



## Proven Architecture

Widely published, industry leading platform built for scalability with better data economics, removing barriers to answering your biggest biological questions

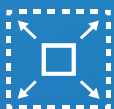

## Flexible Performance

Highly configurable to support the broadest range of methods at any scale – push your research further through cutting edge applications



## Immense Discovery Power

Sequence deeper into the genome, expand into new applications, and run more samples to empower your studies



## Streamlined Operation





Designed to increase lab efficiency with a simplified workflow and seamless user experience



For Research Use Only. Not for use in diagnostic procedures.

illumina®

# NovaSeq System Configurations

				
	S Prime Flow cell	S1 Flow cell	S2 Flow cell	S4 Flow cell
Lanes	2	2	2	4
Output (based on read length)	80–400 Gb	167–500 Gb	417–1250 Gb	2000–3000 Gb
Single Reads (clusters passing filter)	0.8 B	1.6 B	4.1 B	10 B
Run Time	13–38 hours	13–25 hours	16–36 hours	36–44 hours
Max Read Length	2x250	2x150	2x150	2x150
Value†	Smaller batch size Cost effective sequencing for small projects  Pilot new projects QC and optimize library concentrations  Rapid turnaround Faster turnaround time for high-throughput projects‡		Run experiments to scale Enable larger batch sizes and larger cohorts  Attractive per sample economics Competitive pricing per Gb and per M reads  Accessibility to data-rich applications and methods Multi-modal studies, deeper sequencing	


‡ Turnaround time includes cluster generation and sequencing run time

† Compared to Illumina high-throughput portfolio.

# What is a Cluster?



Clusters are a group of DNA strands positioned closely together

Each cluster represents thousands of copies of the same DNA strand in a 1–2 micron spot

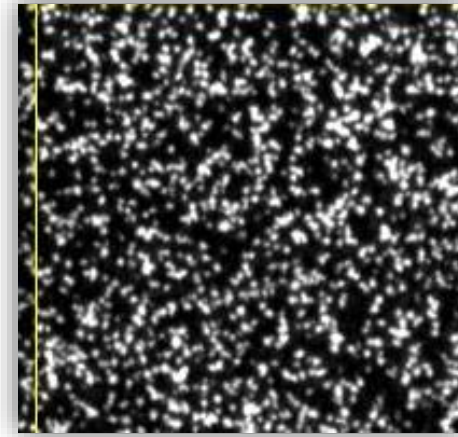

An image of fluorescently labelled clusters on a flow cell

**Single  
DNA  
Library**

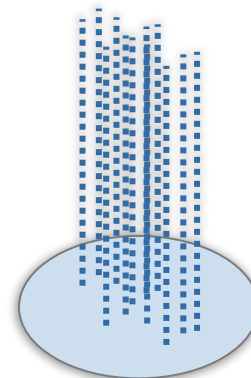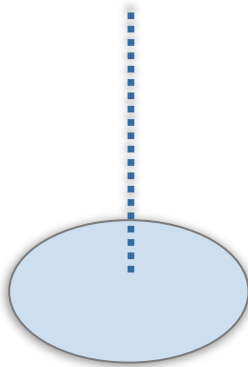



**Amplified  
Clonal  
Cluster**

# What is a flow cell?

Cluster generation occurs on a flow cell

A flow cell is a thick glass slide with channels or lanes

Each lane is coated with a lawn of oligos complementary to library adapters

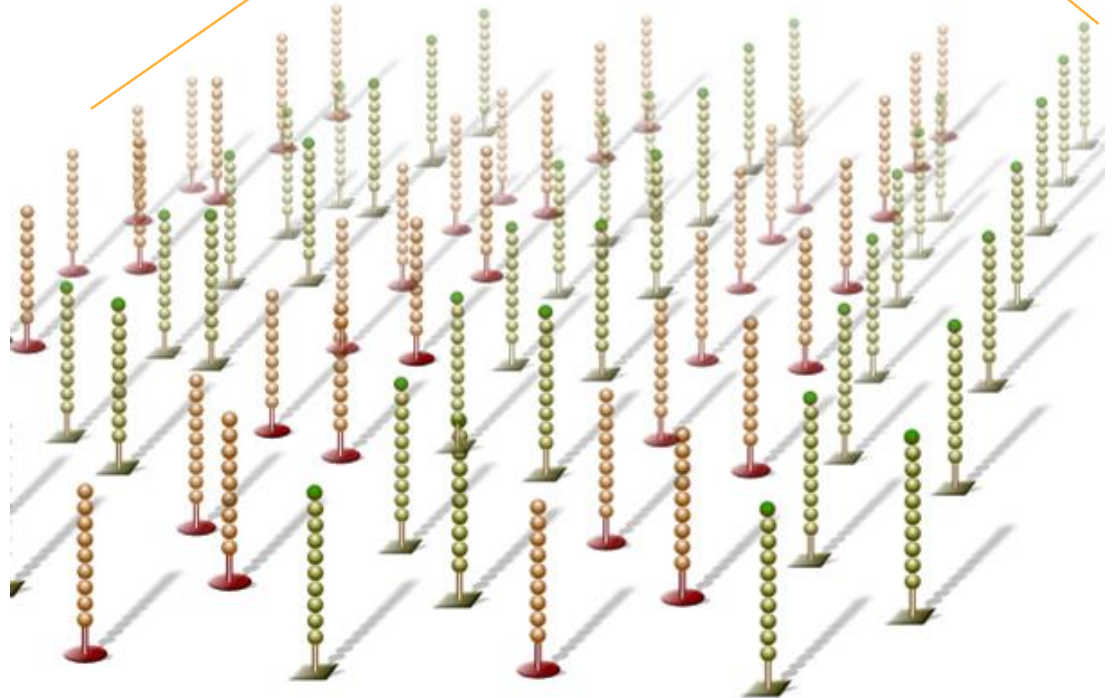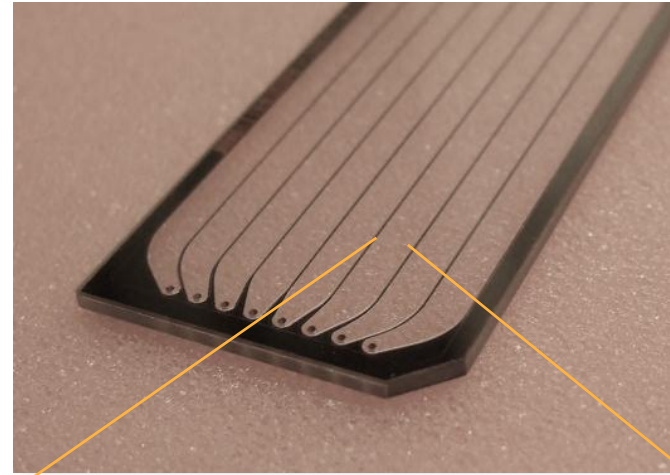


# Flow cell Architecture

## *Random vs Patterned*


### Random Flow Cell

- HiSeq™ 2500, MiSeq™, NextSeq™, MiniSeq™
- Randomly spaced clusters
- Variable Insert Sizes
- Lower Duplication Rates

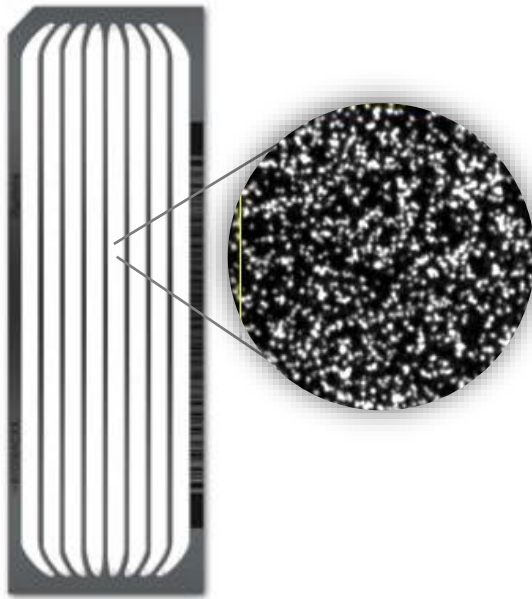


### Patterned Flow Cell

- HiSeq 3K/4K/X, NovaSeq™ 6000, iSeq™ 100
- Defined size and spacing
- Increased Cluster density
- Simplified imaging

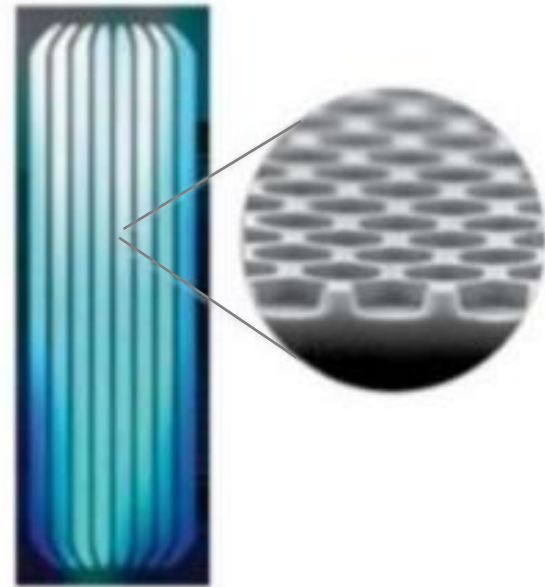


# Patterned flow cells





Complete control of pitch & feature size

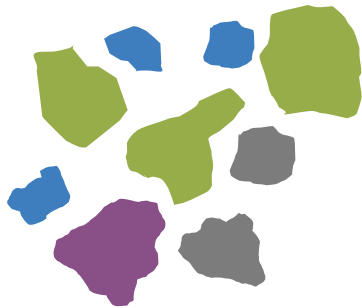

Pitch (Center to Center)

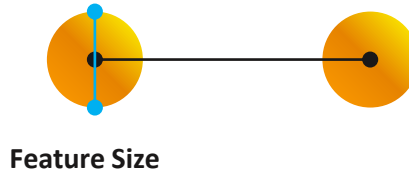

Decrease Pitch

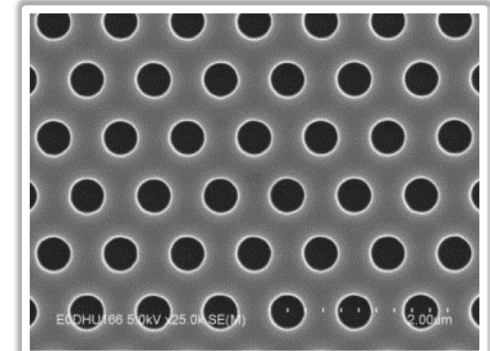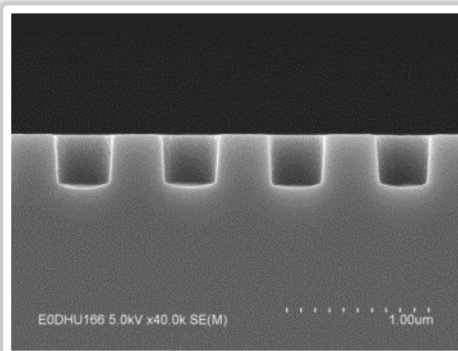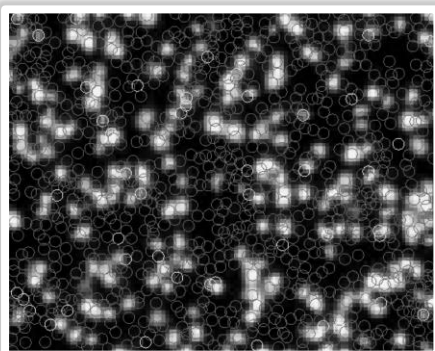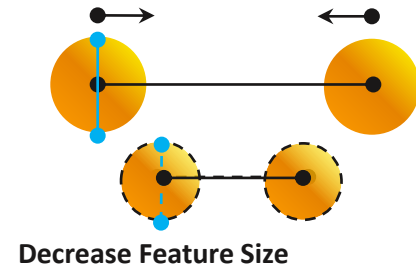

Random

Patterned

Rigid registration reduces time by skipping template generation



# Traditional Cluster Generation





# Hybridize Fragment & Extend


Single-stranded DNA libraries are hybridized to primer lawn

Bound libraries are then extended by polymerases

Surface of flow cell coated with a lawn of oligo pairs

Single-stranded DNA library template molecule

3' extension

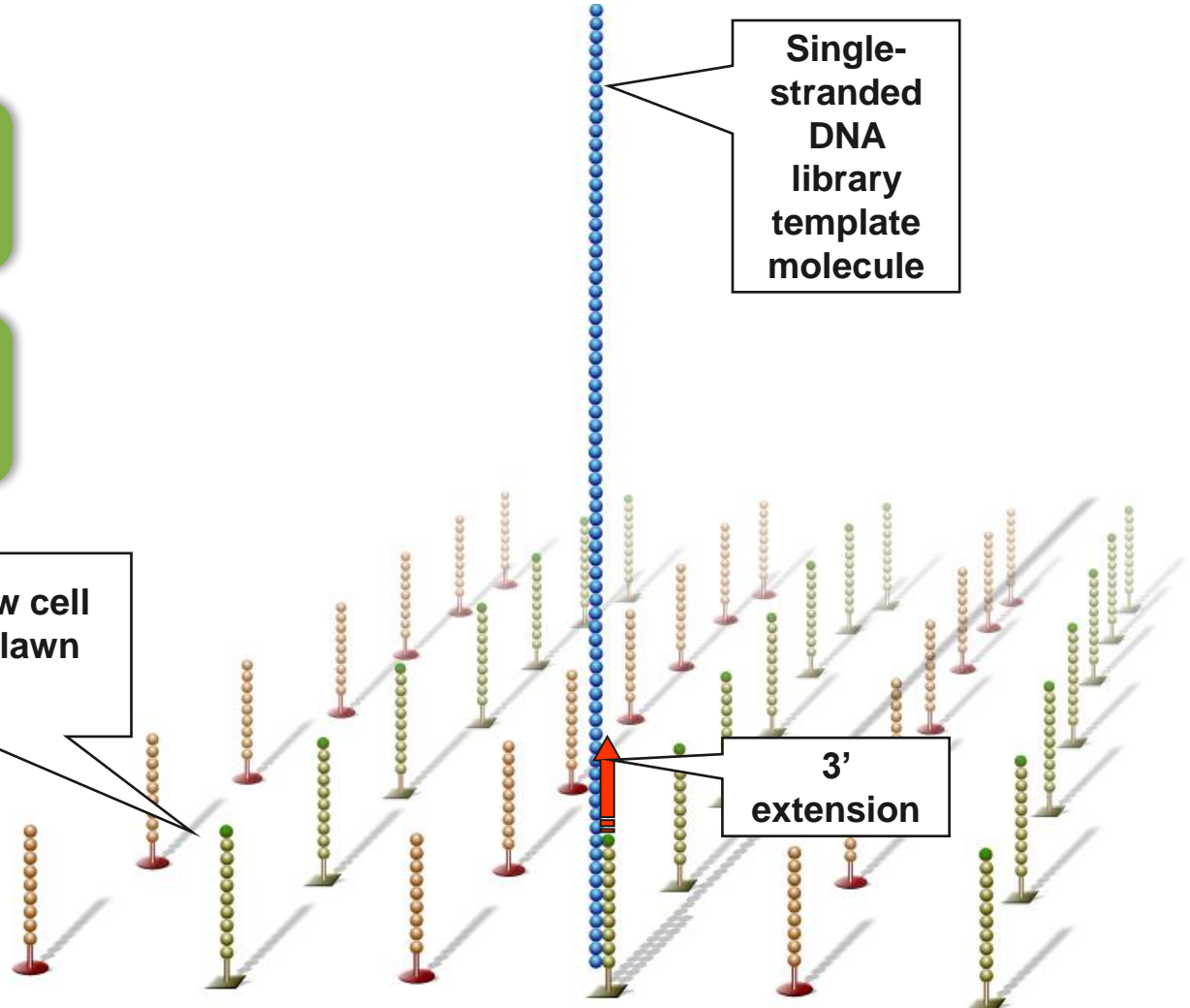


# Denature Double-Stranded DNA


Double-stranded molecule is denatured

Original template washed away


Newly synthesized strand is covalently attached to flow cell surface

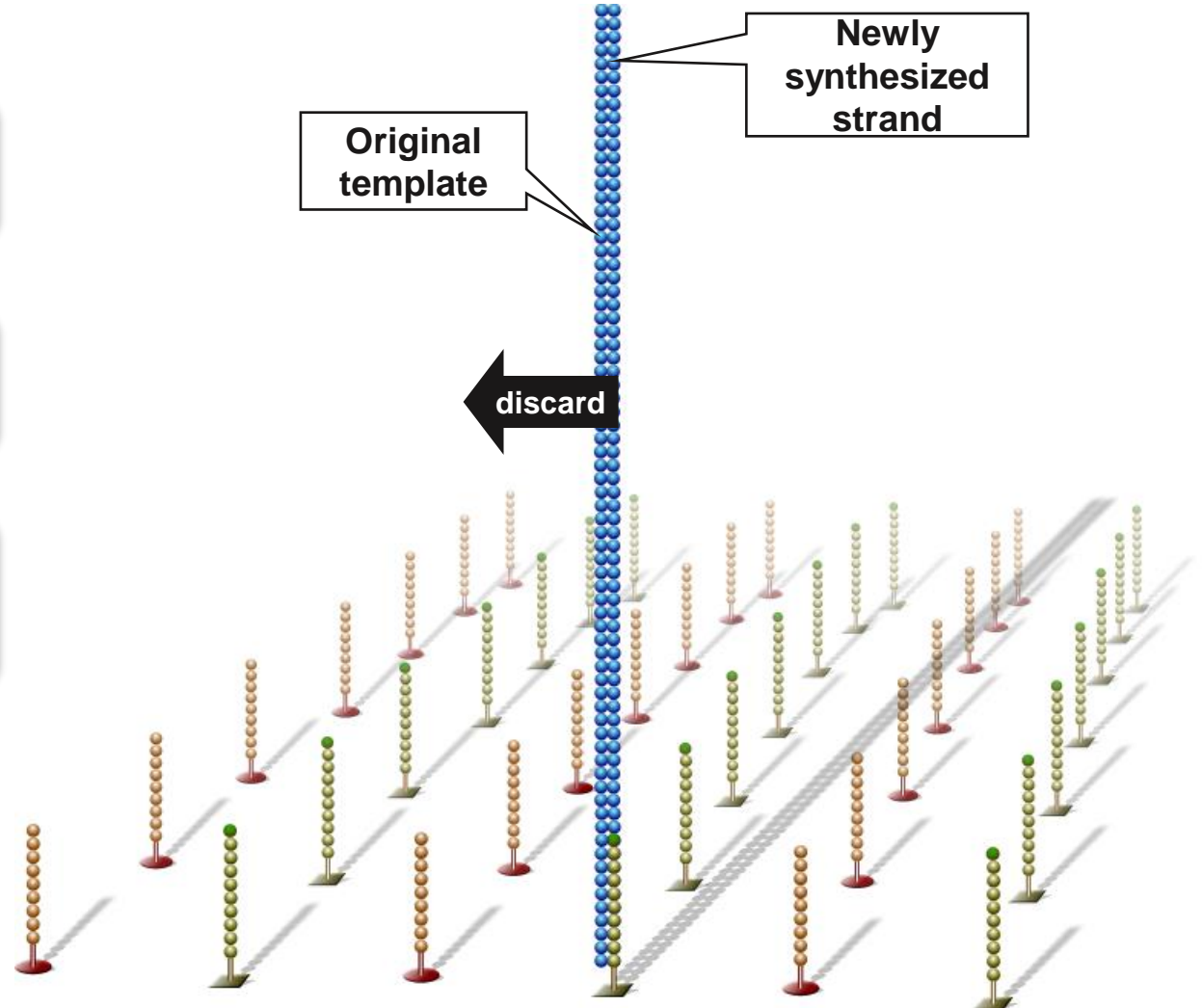

# Single-Stranded DNA

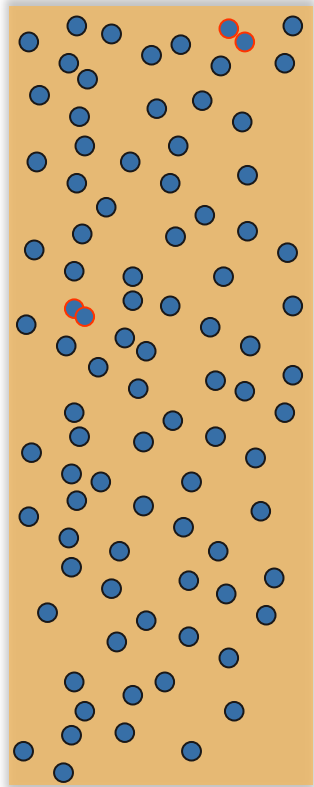


NOTE:  
Single molecules  
bind to flow cell in a  
random pattern

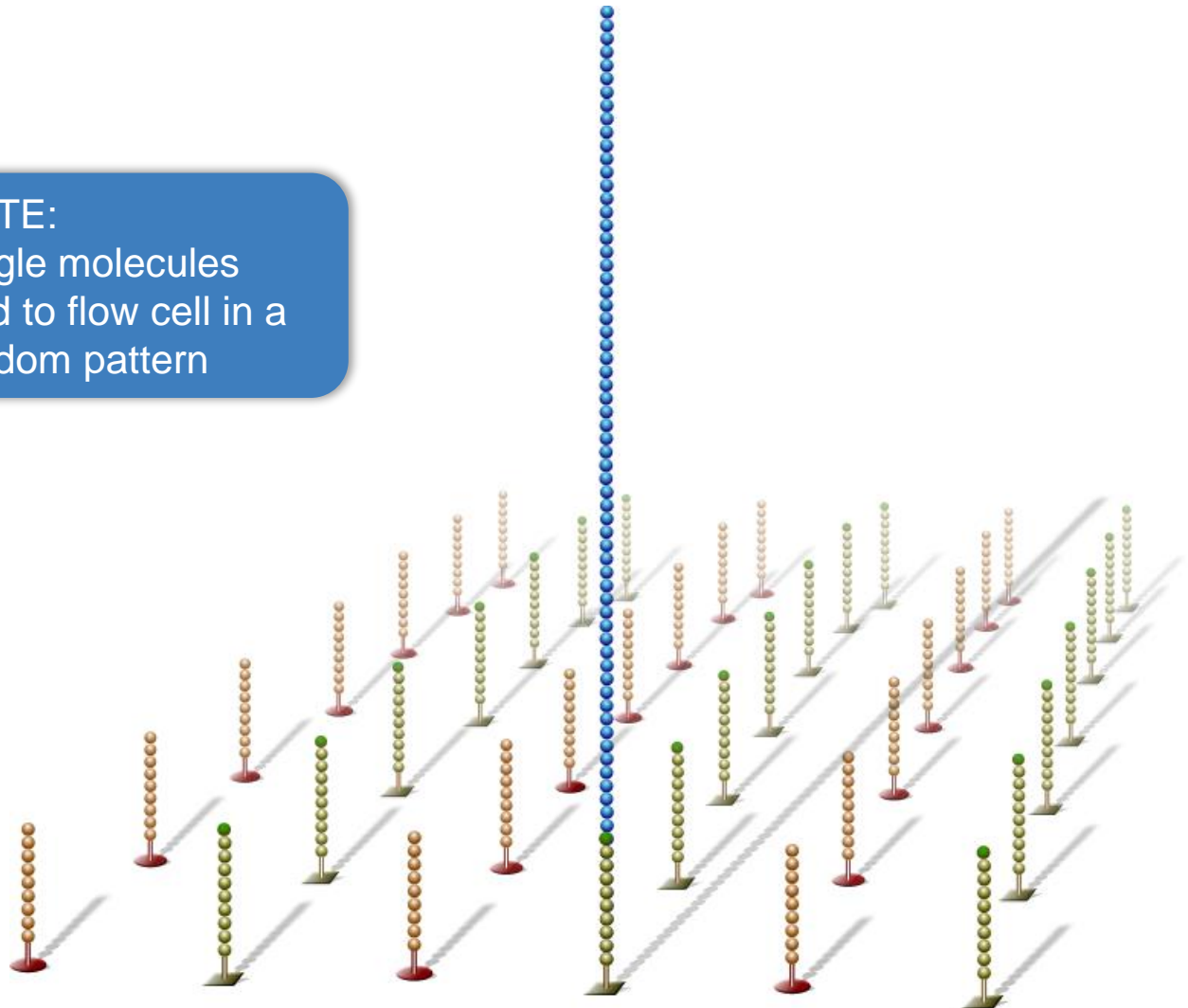

# Bridge Amplification


Single-stranded molecule flips over and forms a bridge by hybridizing to adjacent, complementary primer

Hybridized primer extends by polymerases

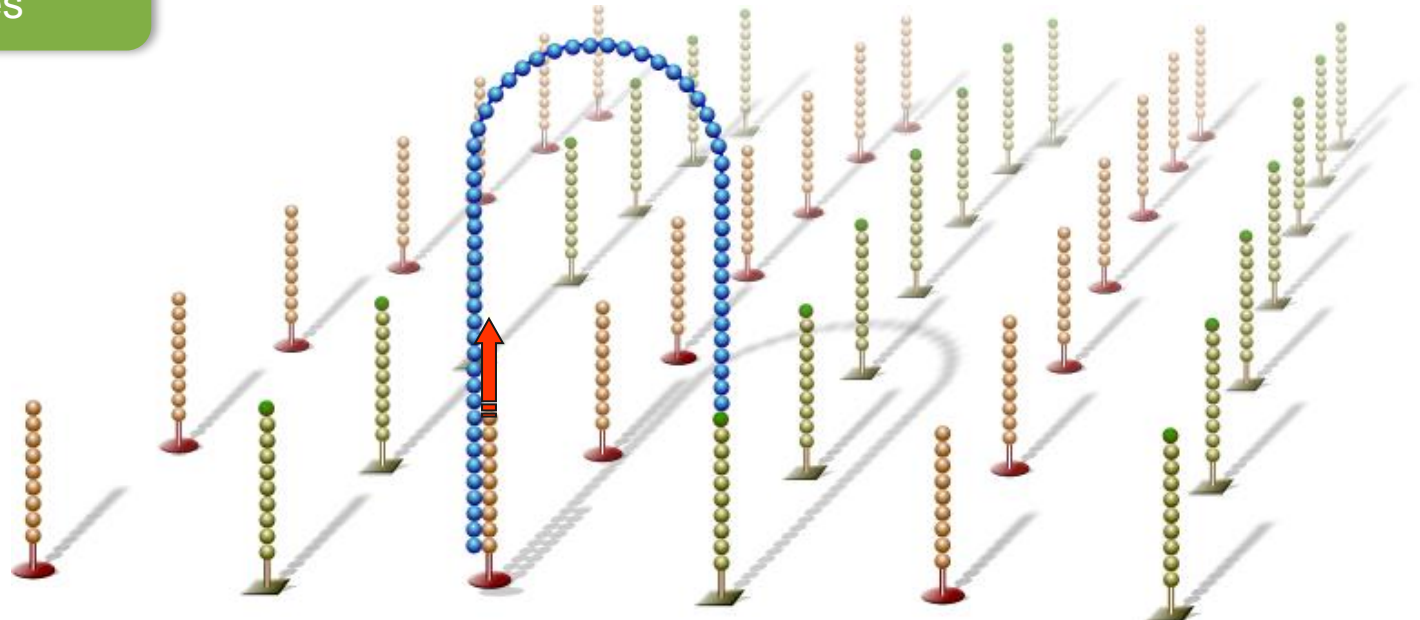

# Bridge Amplification

Double-stranded bridge is formed

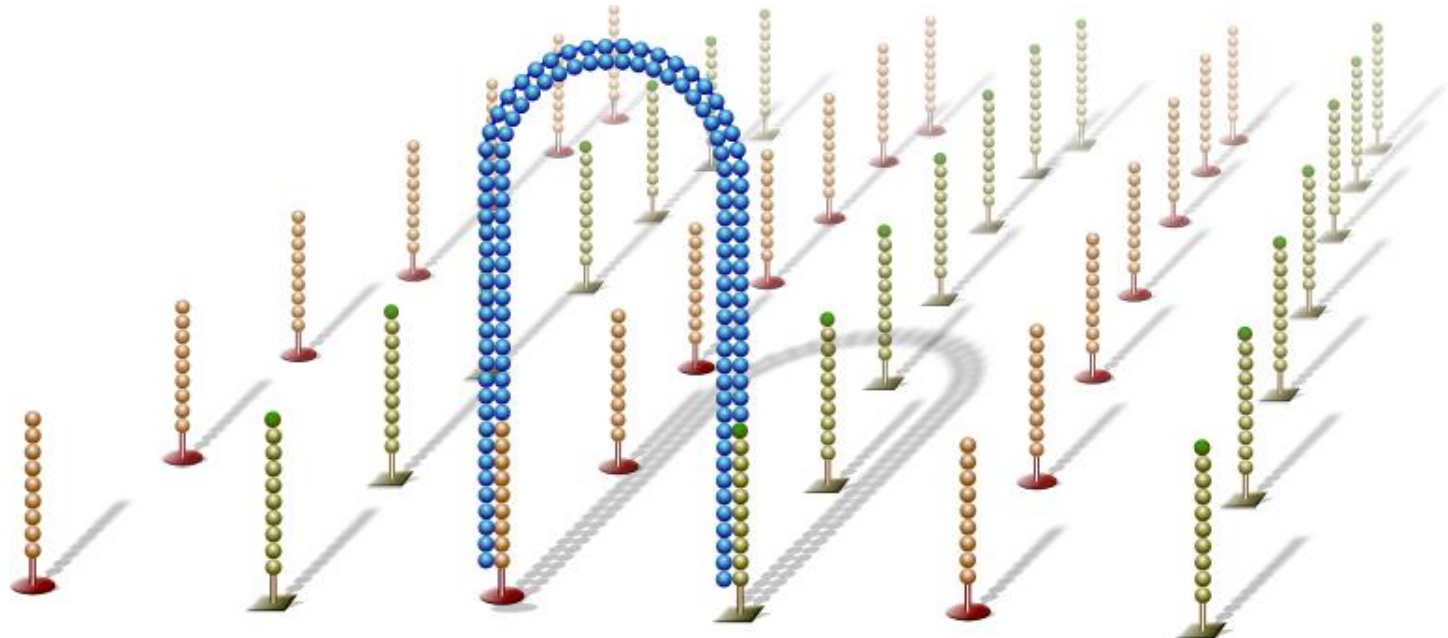


# Denature Double-Stranded Bridge

Double-stranded bridge is  
denatured

## Result:

Two copies of covalently bound  
single-stranded templates

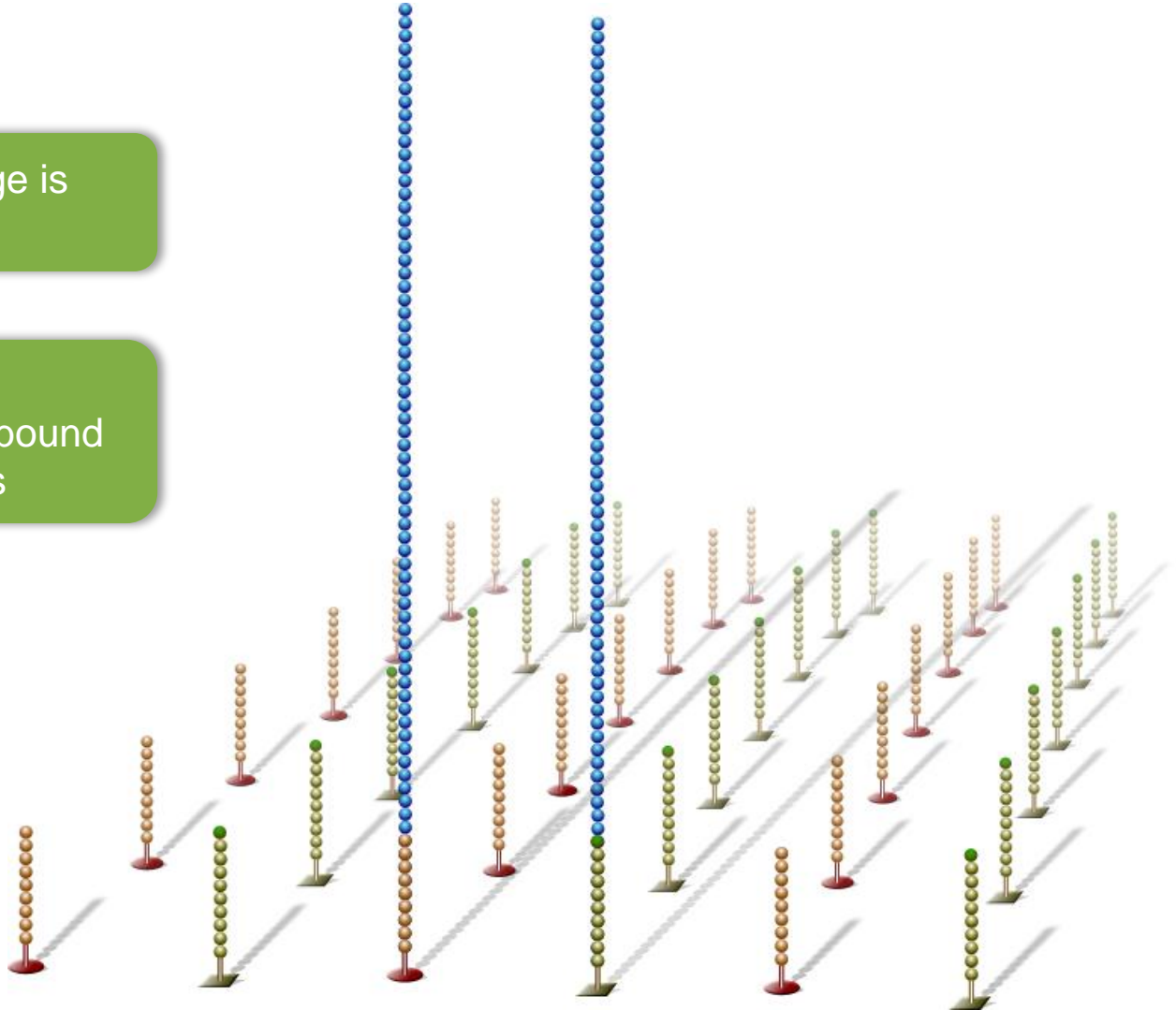




# Bridge Amplification


Single-stranded molecules flip over to hybridize to adjacent primers

Hybridized primer extends by polymerase

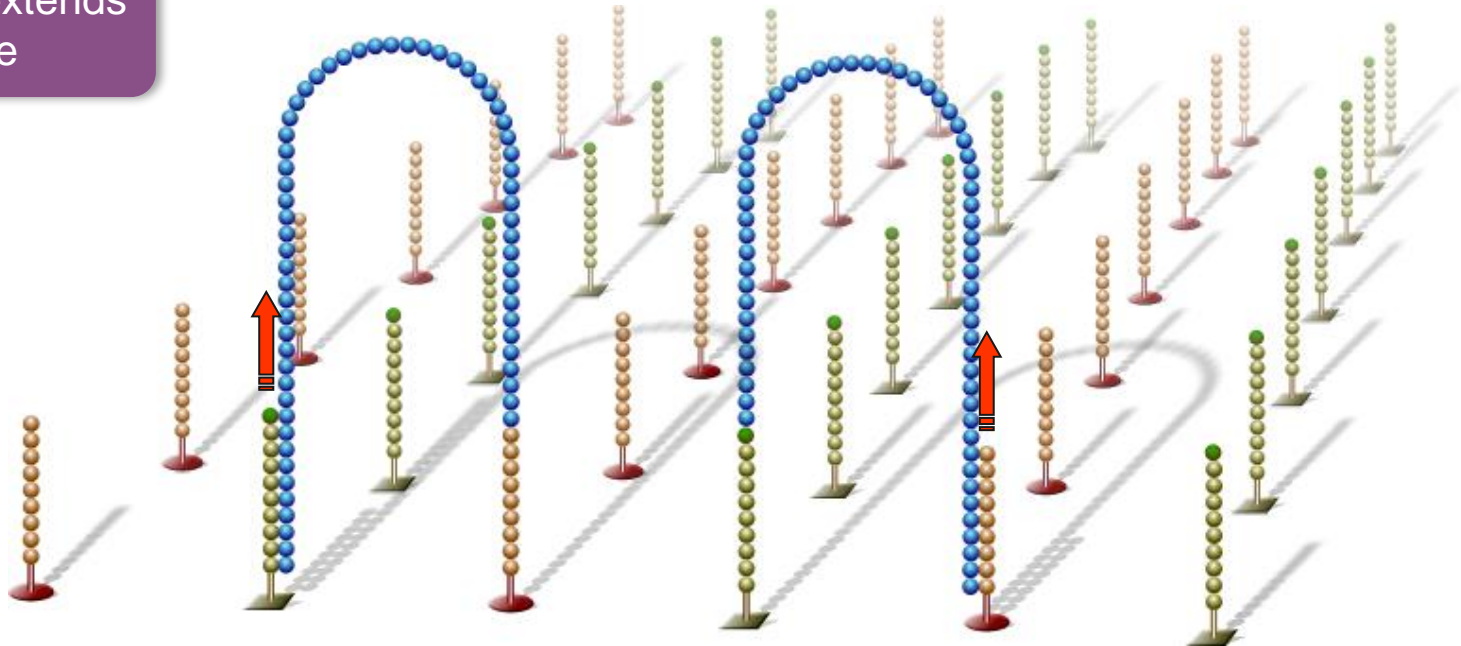

# Bridge Amplification


Bridge amplification cycle is repeated until multiple bridges are formed

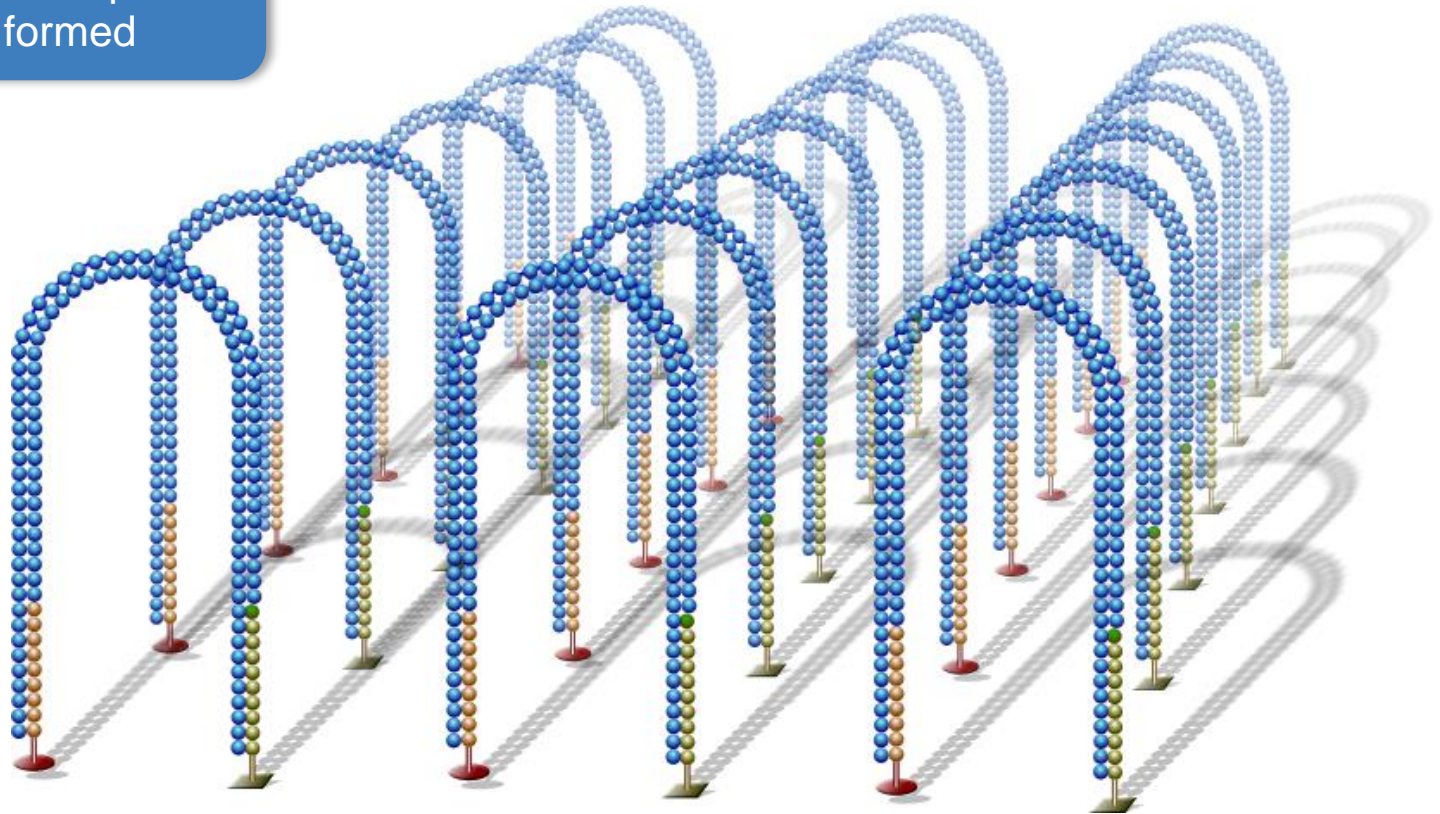




# Linearization

dsDNA bridges are  
denatured

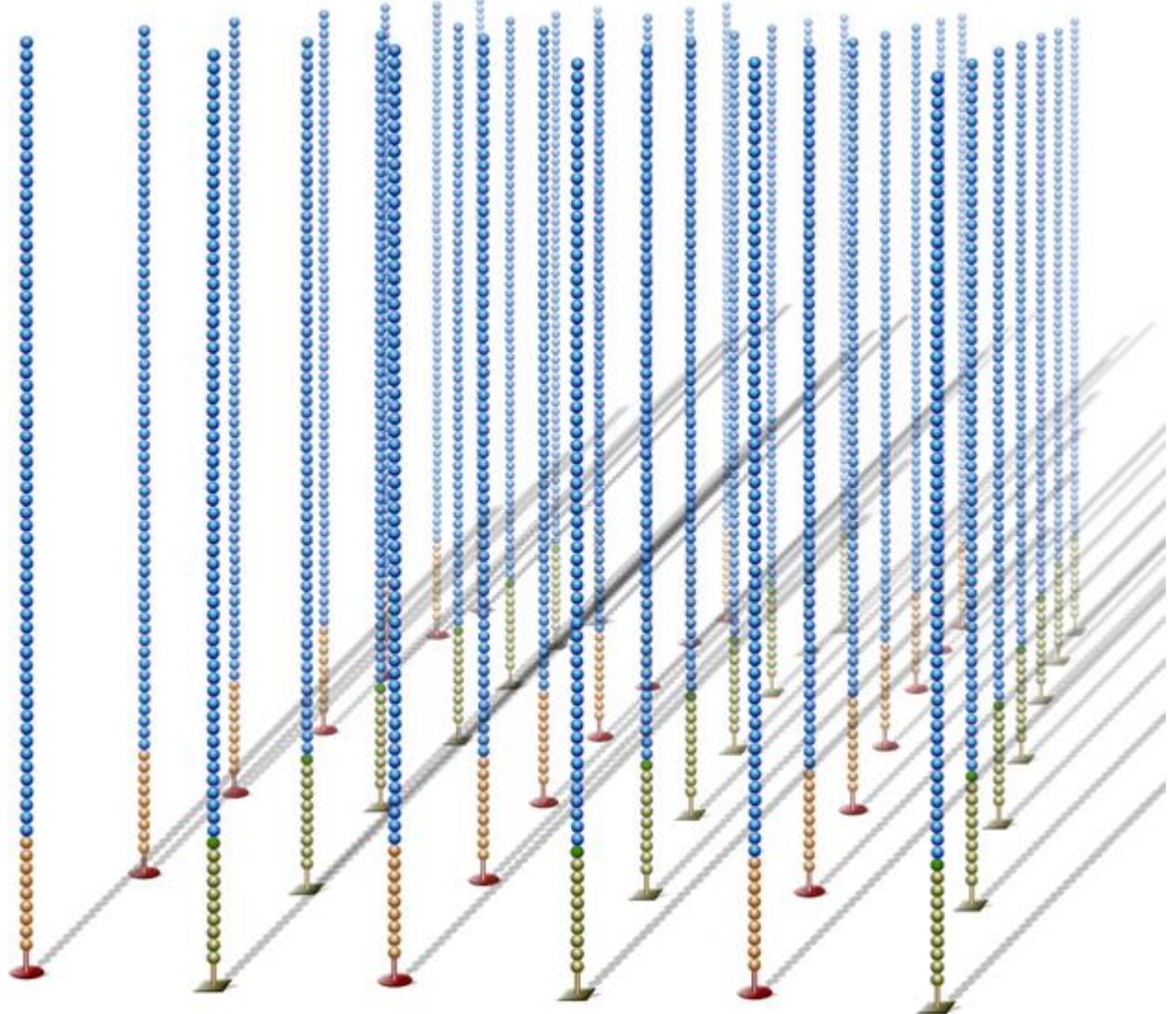

# Reverse Strand Cleavage


Reverse strands are cleaved and washed away, leaving a cluster with forward strands only

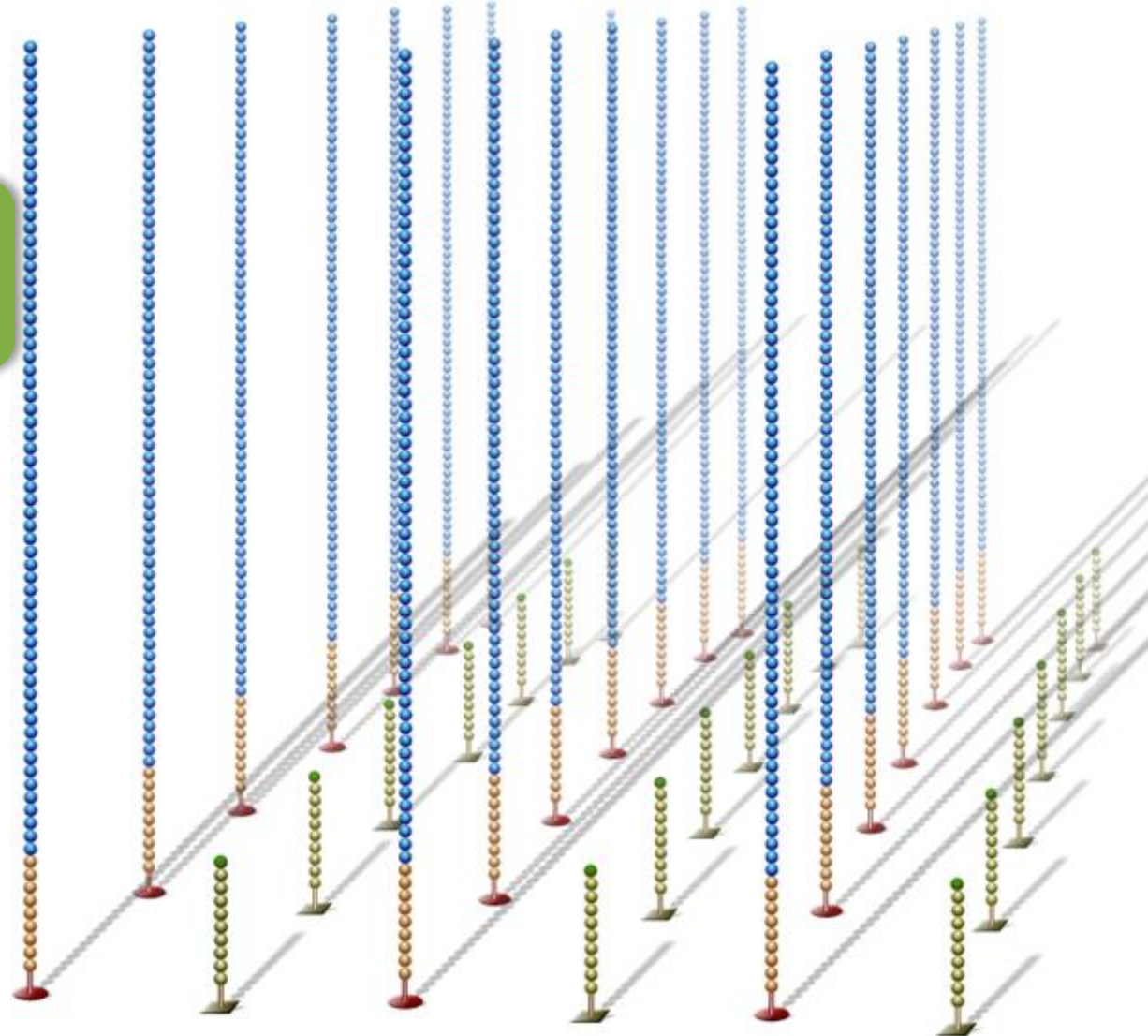

# Read 1 Primer Hybridization


Sequencing primer is hybridized to Read 1 sequencing primer binding site

Sequencing primer

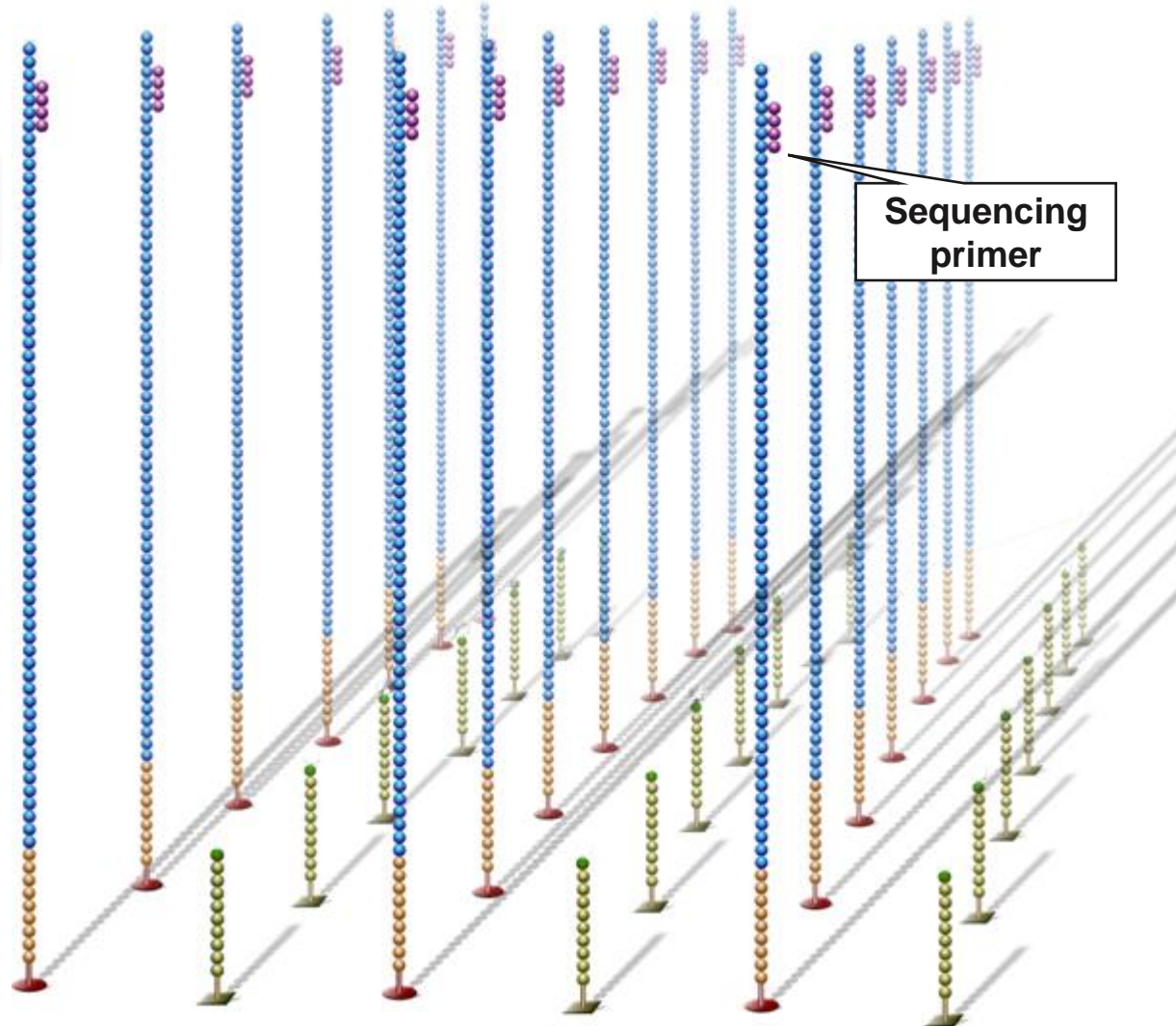

# Illumina Sequencing Workflow


## Sequencing

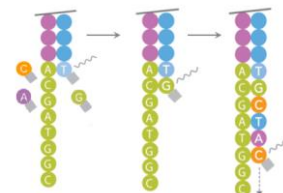[illegible]



# 4-Channel SBS Chemistry:


*HiSeq, MiSeq*

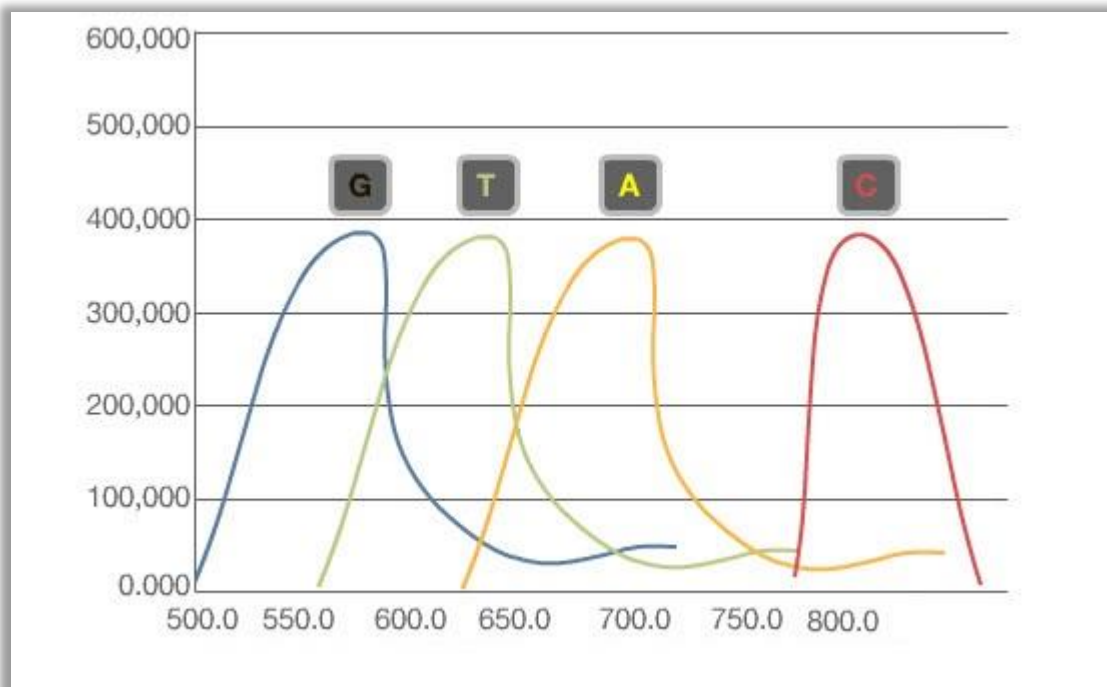

Each of the four DNA bases emits an intensity of a unique wavelength

Collects four images:

- During each cycle, each cluster appears in only one of four images


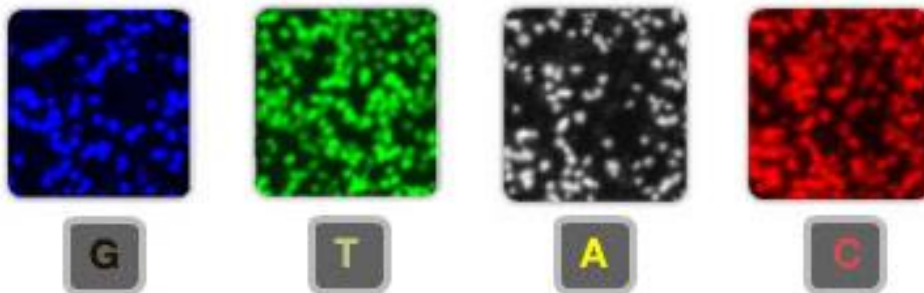


For Research Use Only. Not for use in diagnostic procedures.

# 2-Channel SBS Chemistry:

*NextSeq 550, MiniSeq, NovaSeq 6000*


- 2-channel SBS uses two images:
- Clusters appearing in green only are **T**
- Clusters appearing in red only are **C**
- Clusters appearing in both images are **A**
- Clusters not present/dark are **G**
- After imaging, cluster intensities are plotted and bases called accordingly

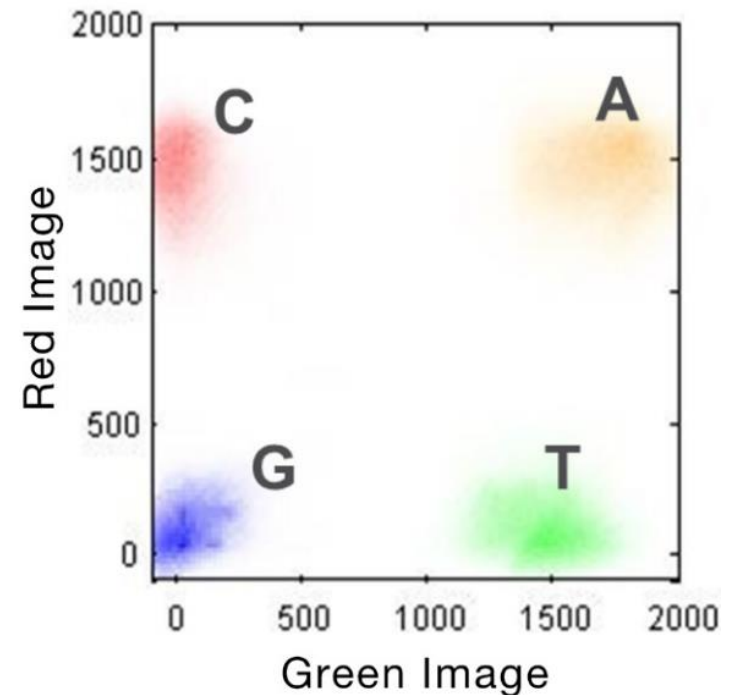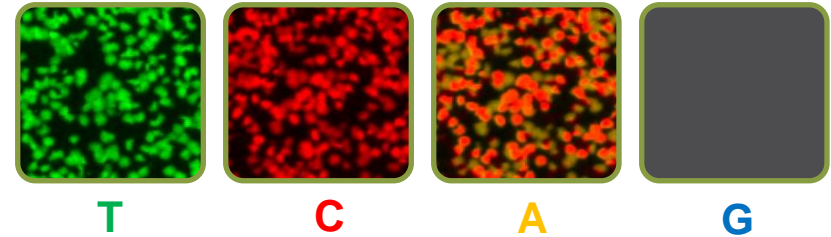

[Illumina Two-Channel SBS Sequencing Technology Technote](#)

# A Closer Look At 2-Dye Chemistry

## *2-channel chemistry*

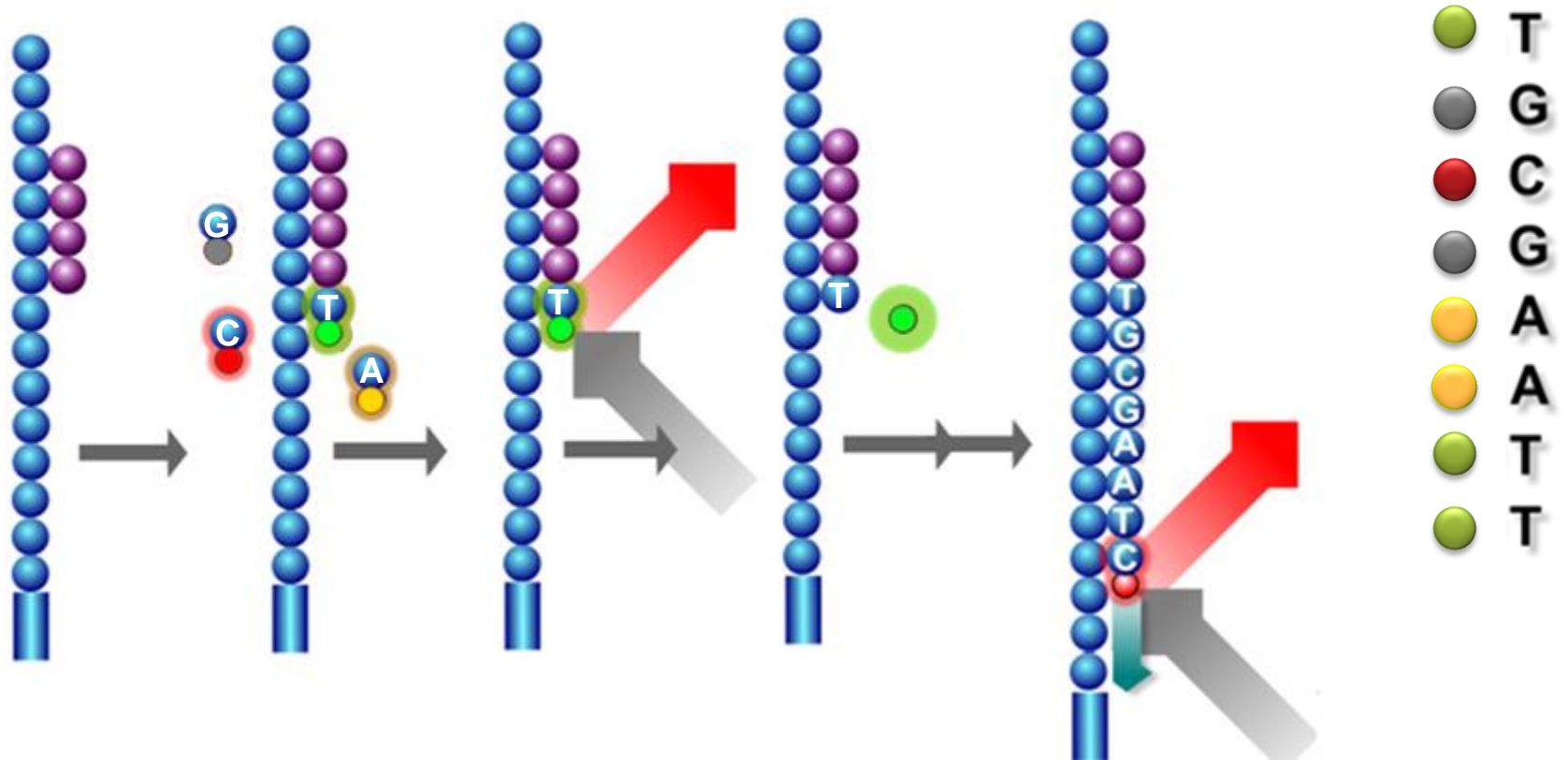









Add 3 FI-NTP's &  
nonFI-NTP + Polymerase








Incorporated FI-NTPs excited  
then imaged (2 images/cycle)








Terminator & fluorescent  
dye cleaved from FI-NTP

X # of bp in read

# Illumina Chemistry Comparison

4-Channel Chemistry				
	 <b>A</b>	 <b>G</b>	 <b>T</b>	 <b>C</b>
Image 1				
Image 2				
Image 3				
Image 4				
Result	<b>A</b>	<b>G</b>	<b>T</b>	<b>C</b>

2-Channel Chemistry				
	 <b>A</b>	<b>G</b>	 <b>T</b>	 <b>C</b>
Image 1				
Image 2				
Result	<b>A</b>	<b>G</b>	<b>T</b>	<b>C</b>

1-Channel Chemistry				
	 <b>A</b>	<b>G</b>	 <b>T</b>	 <b>C</b>
Image 1				
Image 2				
Result	<b>A</b>	<b>G</b>	<b>T</b>	<b>C</b>

----- Intermediate chemistry step

## 4-channel SBS

- Bases are identified using four different fluorescent dyes, one for each base and four images per sequencing cycle

## 2-channel SBS

- Simplified nucleotide detection by using two fluorescent dyes and two images to determine all four base calls

## 1-channel SBS

- Base calling uses one fluorescent dye and two images, with chemistry step in between, to determine all four base calls



# Paired-End Sequencing



# Sequencing with Paired-Ends



**Reference**

This is really the best way to do sequencing

**Single-reads**

This is

...

is really

...

really the

...

the best

...

sequencing


**Paired-reads**

This is (----100 characters-----) sequencing

***Assembly becomes easier***

# Paired-End Sequencing


The read product is stripped off

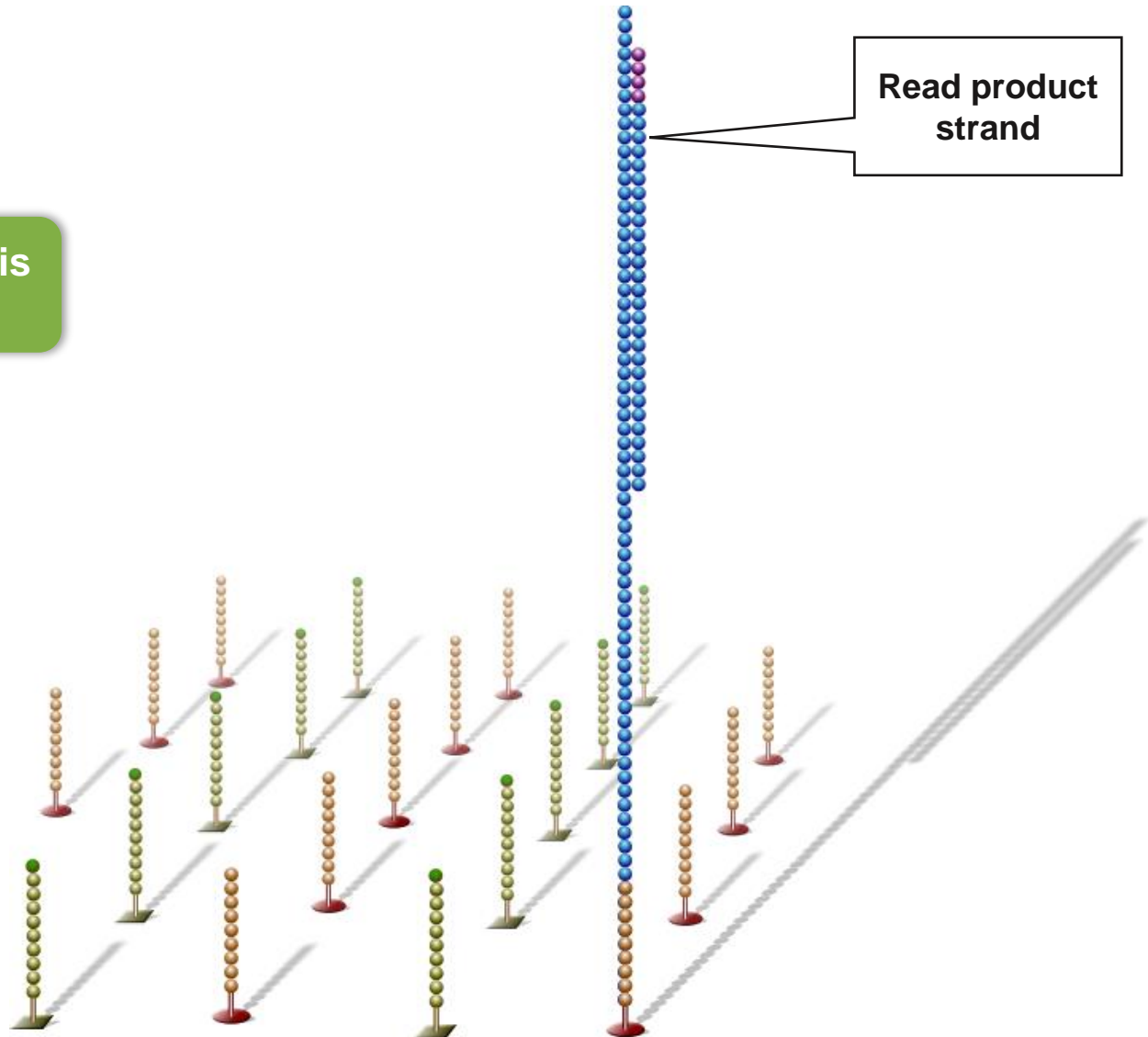


# Paired-End Sequencing

Single-stranded template loops over to form a bridge by hybridizing with a lawn primer

3'-ends of lawn primer are extended

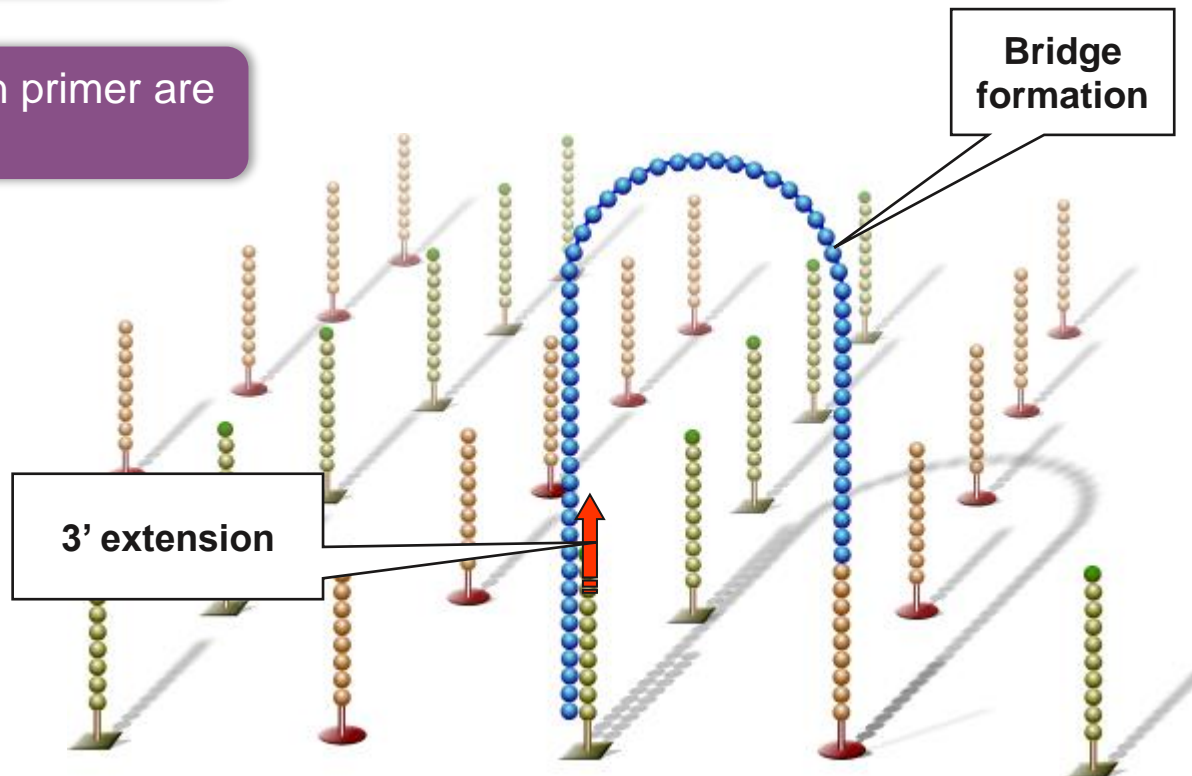


# Paired-End Sequencing

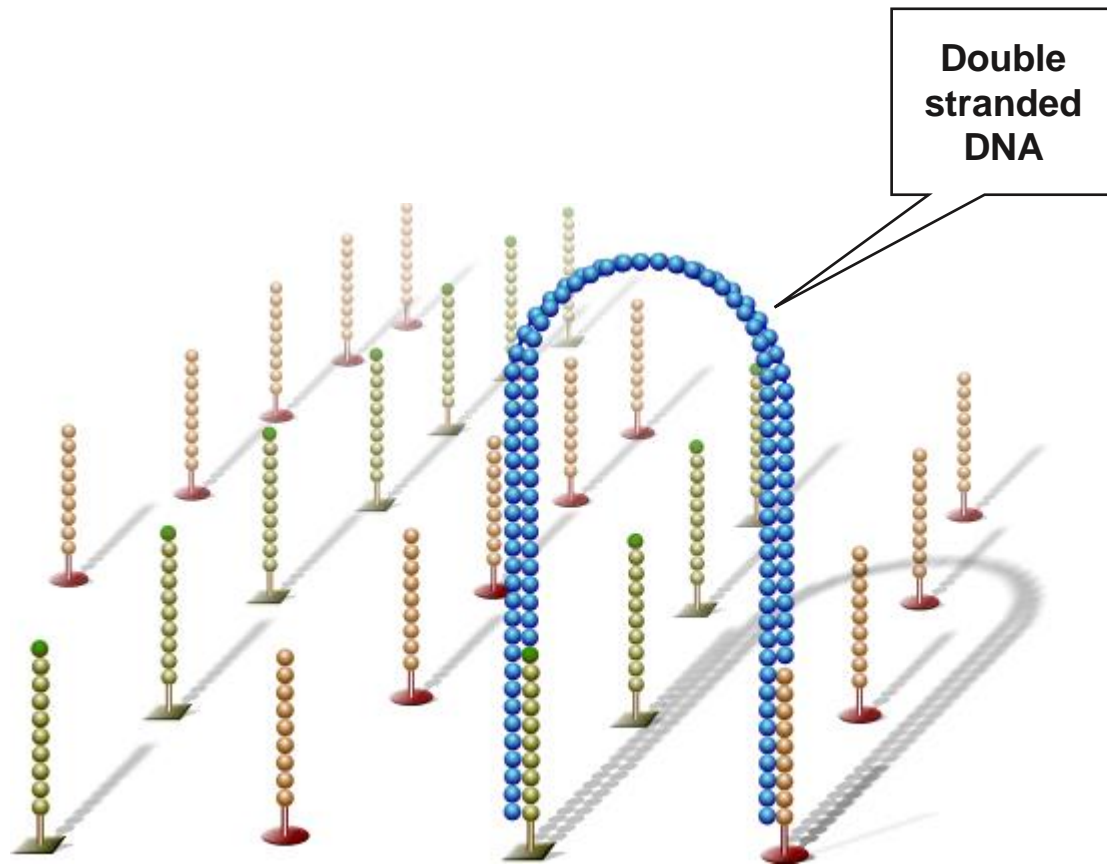


# Paired-End Sequencing

Bridges are linearized and the original forward template is cleaved

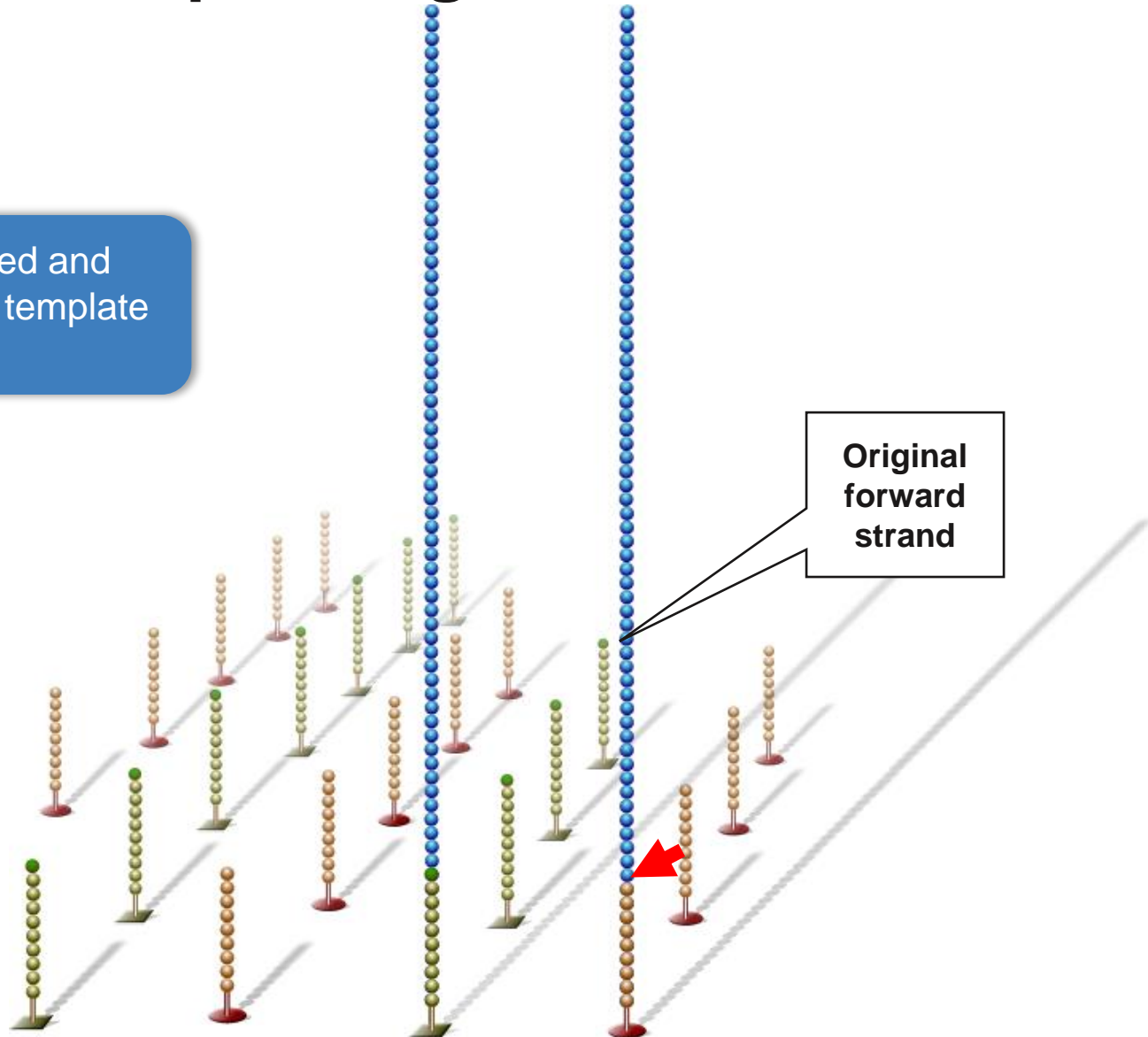


# Paired-End Sequencing

Sequencing primer is hybridized to the Read 2 sequencing primer binding site

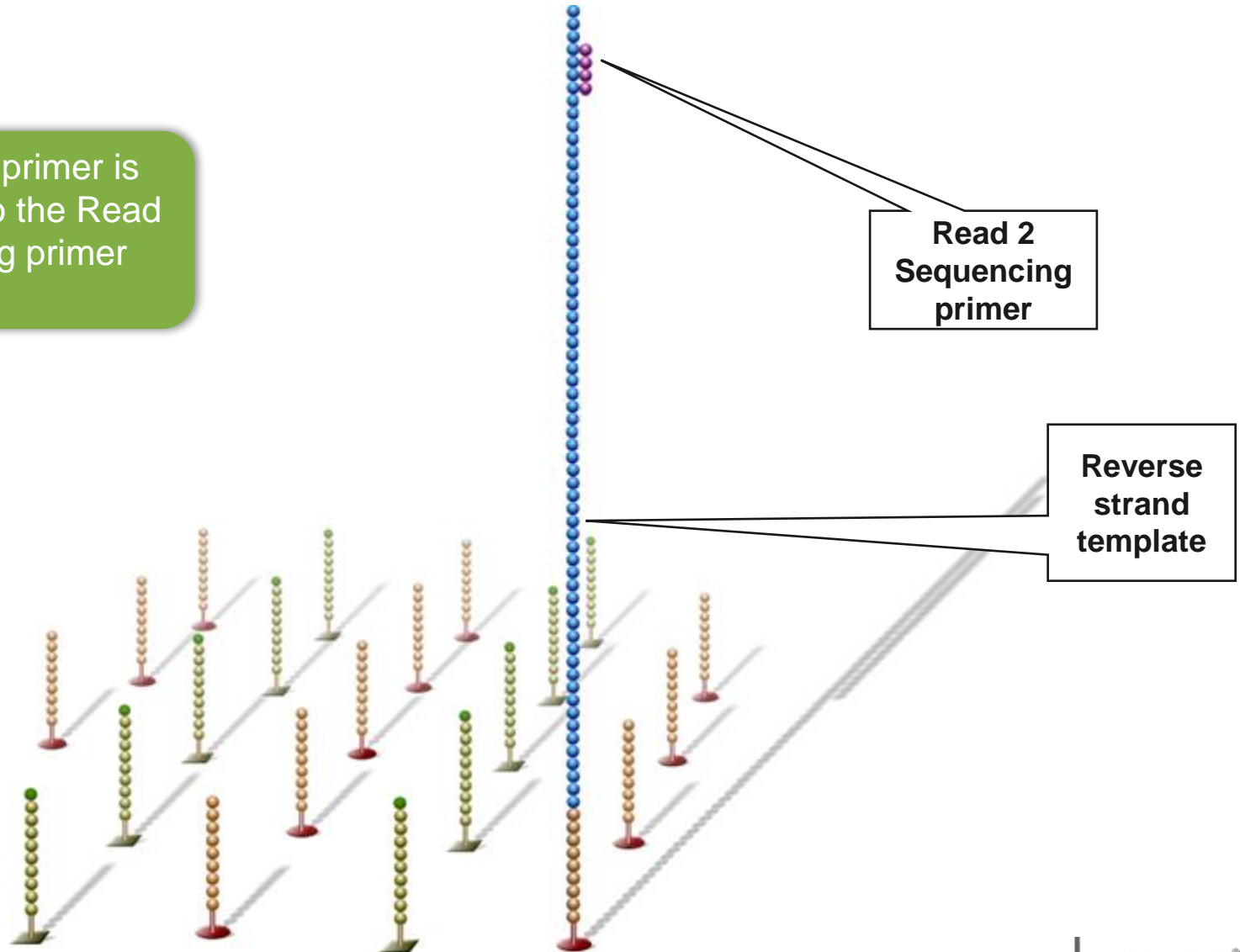


# Sequencing with Index Reads

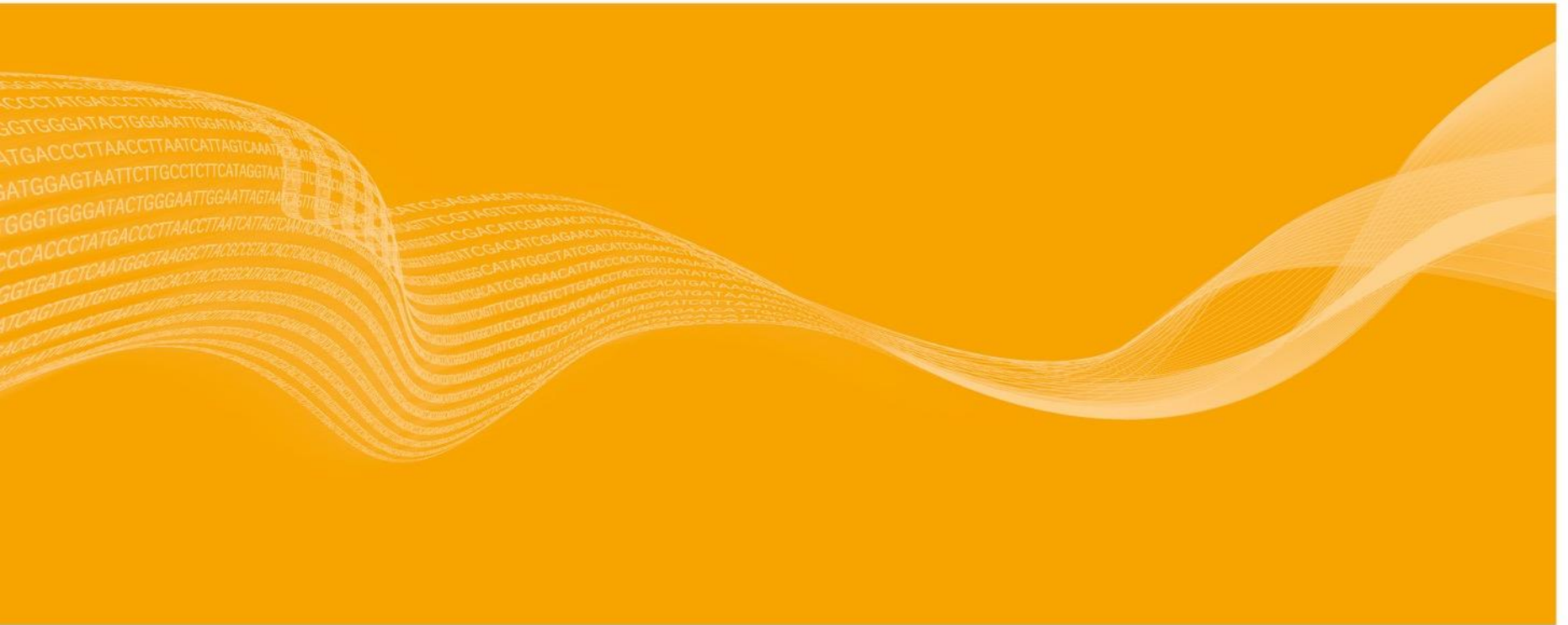




# Single Index Reads

*All Platforms*

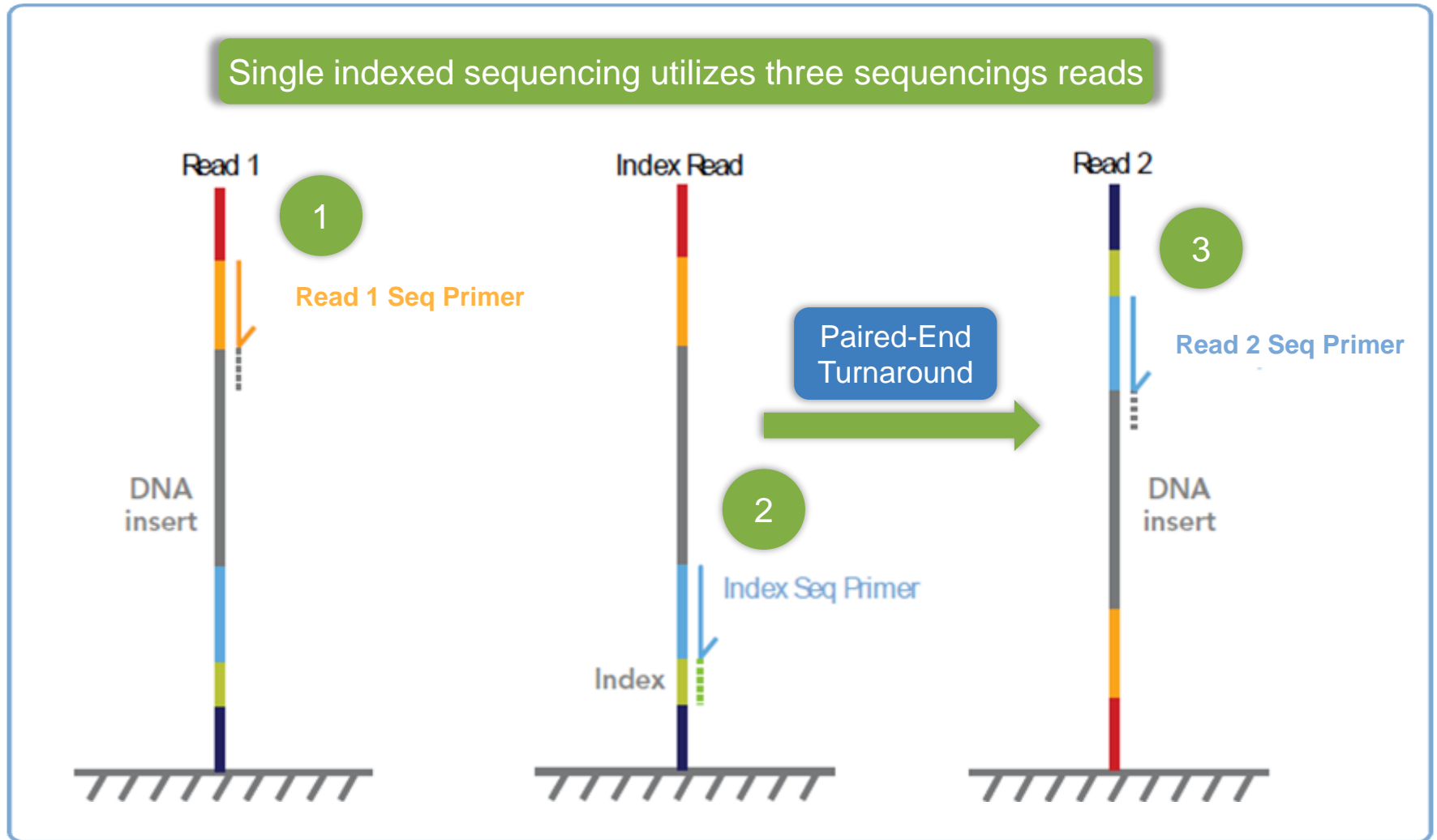

# Dual Index Reads


*MiSeq, HiSeq 2500, NovaSeq 6000*

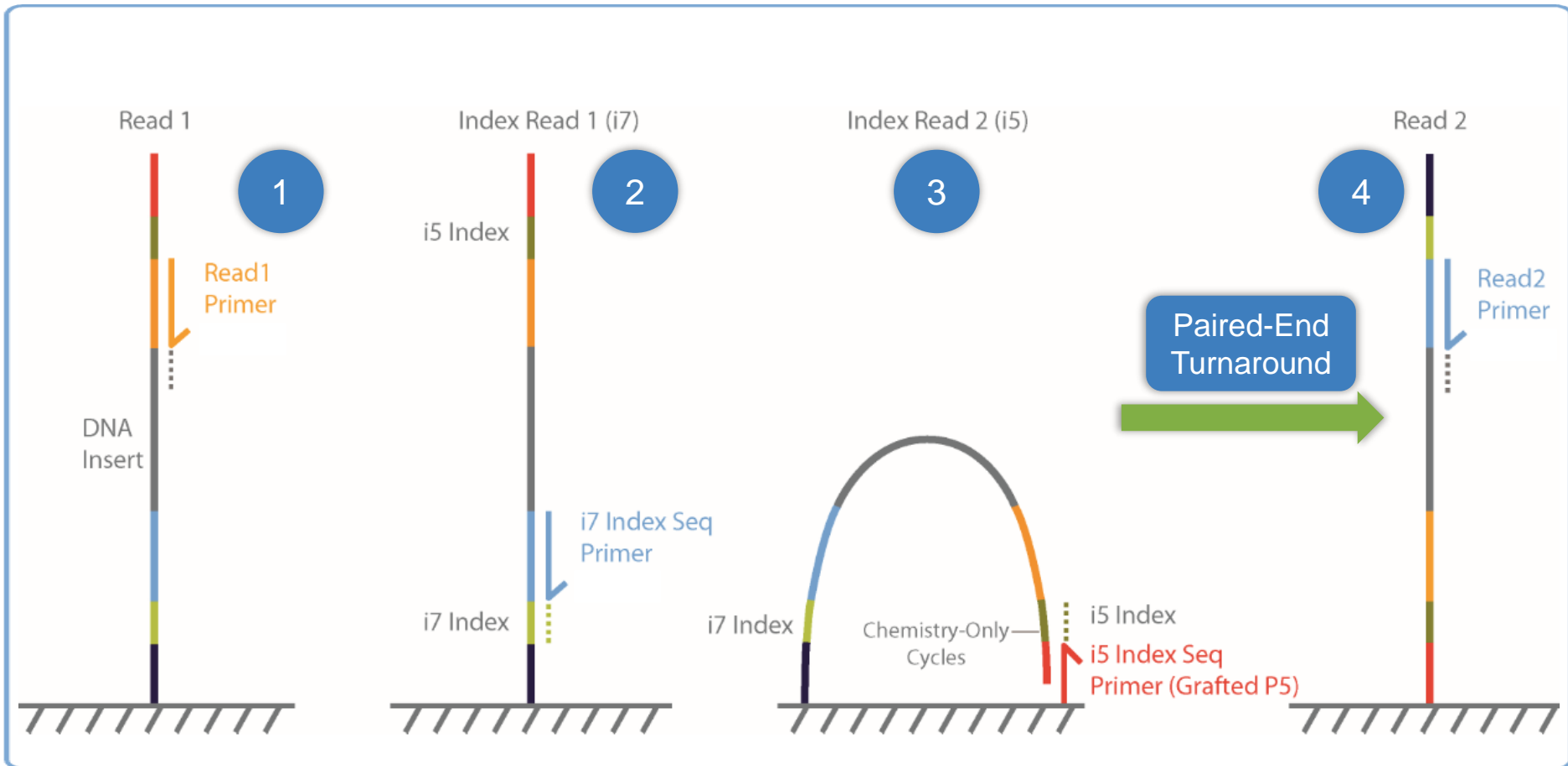

Dual indexed sequencing utilizes four sequencing reads

# Dual Index Reads

*iSeq 100, MiniSeq, NextSeq, HiSeq 3000/4000*

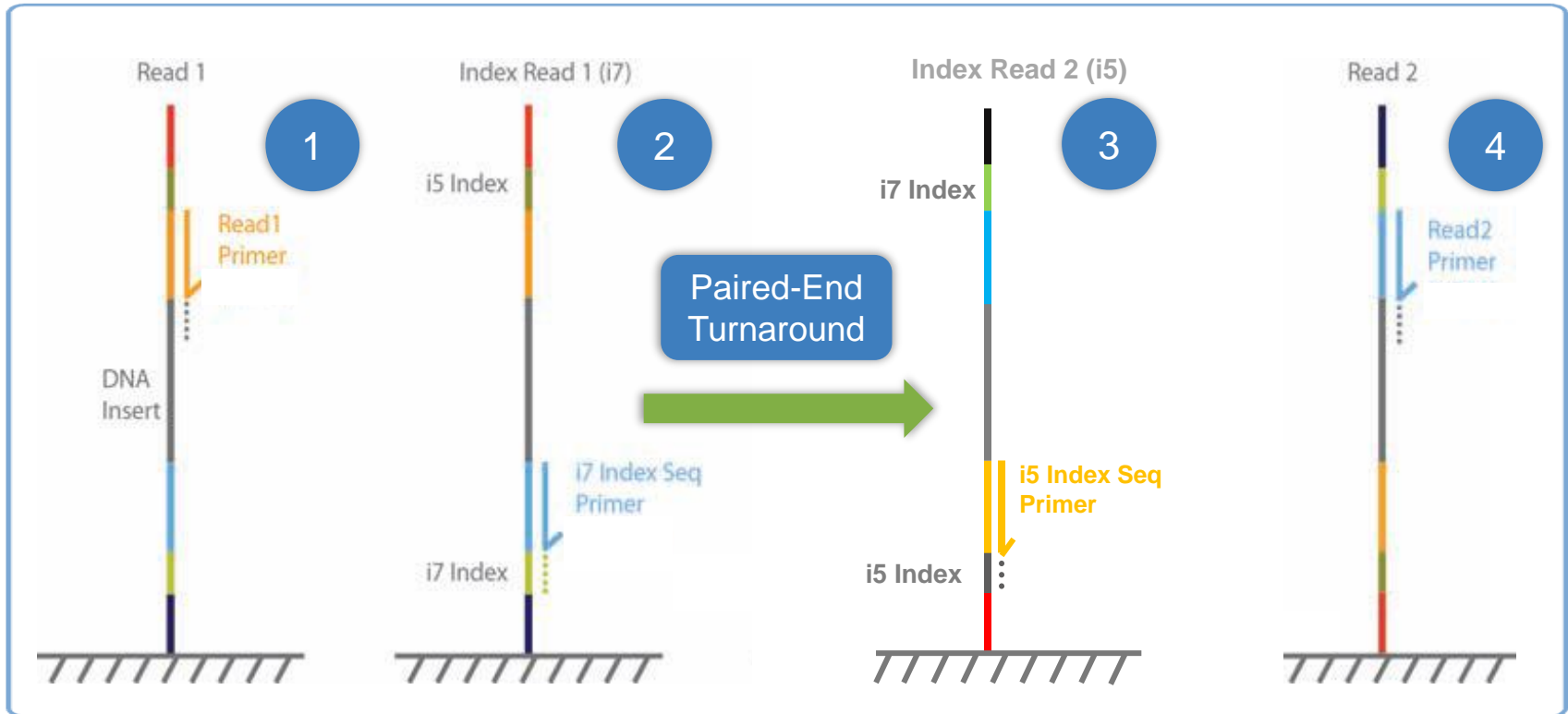


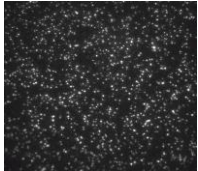




# Illumina Sequencing Workflow

## Sequencing

## Data Analysis

STAAGGCTAGGTTTCATGCTA  
STAAGGCTAGGTTTCATGCTA  
STAAGGCTAGGTTTCATGCTA  
STAAGGCTAGGTTTCATGCTA  
TAAGGCTAGGTTTCATGCTA  
AAGGCTAGGTTTCATGCTA  
ST AGGCTAGGTTTCATGCTA  
STA GCTAGGTTTCATGCTA  
STAA CTAGGTTTCATGCTA

# Analysis Overview

Analysis Type	Software	Outputs
Control Software		 Images, Intensities and Base Calls
Analysis Software		 Alignments, Variant Detection
Visualization Software		 Annotation, Filtering, Reports

# Questions?

Ryan Gentry M.S.

[rgentry@illumina.com](mailto:rgentry@illumina.com)

813-404-1082

